



EXPERIENCED  
HUMAN  
SOLUTION-  
FOCUSED



STIEGELMEYER  
GROUP

# Catalogue



Care beds and accessories  
for private customers

---

# Contents

---

## General information

- 3 Burmeier, the company
- 4 Safety standards for care beds
- 5 The 24-volt drive system

## Care beds, bed insert frames and other products

- |                               |                          |
|-------------------------------|--------------------------|
| 6 Regia                       | 34 Lippe IV              |
| 10 Regia partner              | 38 Lippe 120/140         |
| 14 Regia gardo                | 42 Relax                 |
| 18 Inovia                     | 46 Gigant                |
| 22 Westfalia-Klassik          | 50 Mattresses            |
| 26 Dali with wooden surrounds | 52 Bedside cabinets      |
| 30 Dali economic              | 55 The Stieglmeyer Group |



## The company

# Safe living at home

---

More support for caring for people at home – this goal underpins our work at Burmeier. The company, with its roots in Hiddenhausen in East Westphalia, is the market leader in the home care sector in Germany. The technically innovative and durable products enjoy a worldwide reputation.

In its current form, Burmeier has been another successful subsidiary of the Stieglmeyer Group since 1995. Product development always focuses on the concept of making

life at home as pleasant as possible for people in need of care while also facilitating the work of care-giving relatives.

An attractive design, modern safety side systems and prompt availability of the products are key arguments for retailers and private customers. Discover on the pages ahead how Burmeier care beds can decisively boost your quality of life in a home care setting.





## Safety standards

All beds are designed for maximum safety as standard.

### This includes:

- Primary fuse, electronic overload cut-off device and thermal fuse against overheating
- Motors run on safety extra-low voltage
- In-company quality standards exceed industry standard requirements
- Anti-kink strain relief on bed frame and on motor housing
- At least IPX4 splash-proof

The following standards and directives are applied

- Regulation (EU) 2017/745 (Medical Device Regulation)
- DIN EN 62353: Medical electrical equipment – Recurrent test and test after repair of medical electrical equipment
- DIN EN 12182: Technical aids for disabled persons
- DIN 32977: Adapted designing for disabled persons
- DIN EN 60601-1: Medical electrical equipment – General requirements for basic safety and essential performance
- DIN EN 60601-1-2: Medical electrical equipment – Collateral standard: Electromagnetic compatibility
- DIN EN 60601-1-6: Medical electrical equipment – Collateral standard: Suitability for intended use
- DIN EN 60601-2-52: Medical electrical equipment – Particular requirements for basic safety and essential performance of medical beds
- DIN EN ISO 14971: Medical devices – Application of risk management to medical devices

## The 24-volt drive system

Reliable safety, low costs and trouble-free operation: The original 24-volt system from Burmeier allows you to reap the benefits of your care bed even more effectively.

With its drive system, all live 230-volt components are located outside the bed. Bed users enjoy the added safety of the 24-volt safety extra-low voltage, while bed operators benefit from minimised costs: Electricity costs are appreciable lower than with conventional systems.


Even greater costs and effort are saved due to the significant reduction in leakage current measurements. Commercially used care beds normally require to be leakage current tested every 2 years. The use of the 24-volt drive system prolongs this interval so significantly that no further measurements are required within the lifetime of the bed and beyond. The length of this lifetime is given in the operating instructions.

Leakage current testing is not even required during initial commissioning since the manufacturer has already tested the drive. This means that electrically trained persons do not have to be present when the bed is set up.

### Further advantages at a glance:

- Certified by TÜV Rheinland
- Ergonomic handset with selective locking function
- The care bed has no current-carrying components except when the handset is in use
- No electrosmog at the bedside
- Also available in conjunction with Linak double drives and an OpenBus drive system

## Strong advantages

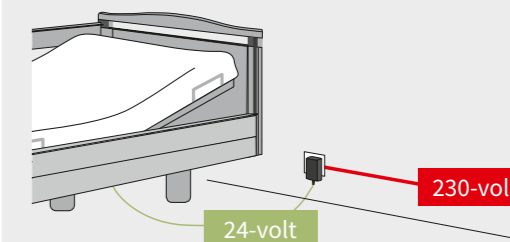

- Safety extra-low voltage starting from the transformer at the power socket
- No 230-volt components on the bed
- No leakage current measurements are necessary after assembling the bed or beyond during the expected service life of the bed
- A test certificate can be downloaded for every bed



*Regia*

# Regia

The care bed that offers many options



### Large height adjustment range

With its large height adjustment range of 26 to 80 cm, the Regia satisfies all the requirements of the care sector. The lowest position is extremely low and plays an important role in preventing fall injuries. With its 3-stop function, the bed automatically stops at an intermediate position of around 40 cm – the ideal height for conveniently getting in and out of bed. Raised to its highest position, the Regia enables caregivers to work ergonomically at hip level.

### Telescopic safety sides

The telescopic safety sides (TSG) supplied as standard are easy to raise in 2 stages. If they are only raised at the head end, this leaves ample room for getting in and out of bed. In this way, Burmeier follows the Werdenfelser approach, which consistently avoids the use of liberty deprivation measures (DoLS).

### Easy Switch system

The modular Easy Switch system makes it possible to attach and remove the safety side elements, the headboard and footboard and the side panels without tools. This enables the bed to be adapted to the changing needs of the resident at any time. Attractive fabric and faux leather upholstery is also available on request for the panels and the headboard and footboard. There is also a quilted option to choose from.

### Safe working load of 225 kg

The Regia also offers enhanced comfort for heavier residents.

### Features

Mattress base: 90 x 200 cm or 100 x 200 cm, in 4 sections, with easy-clean metal slats; mattress base height adjustment from 26 – 80 cm, all mattress base elements electrically adjustable, reverse-Trendelenburg position up to 15° with automatic auto contour setting; mattress base with automatic backrest retraction; variety of headboards and footboards and decors to choose from



▲ Telescopic safety side with Easy Switch system



▲ Reverse-Trendelenburg position



▲ Upholstery with fabric or faux leather\*



▲ Bed extension and mattress extension piece\*



▲ Foam leather cover\*



▲ Linen holder\*



▲ Reading lamp\*



▲ Under bed light\*



▲ Comfort mattress base\*

\*optional accessories



*Regia partner*

# Regia partner

## The double bed for home care



### New flexibility for couples

Accommodating nursing care needs in older age can be a particular challenge for couples. Often, one partner moves into a care bed while the other remains alone in the double bed. The intimacy of the relationship suffers as a result. Our Regia partner offers a solution to this problem: being a proper care bed in itself while also serving as a double bed, it provides couples with the support they need.

### Separates into two single beds

The Regia partner can be separated into two single beds in just a few easy steps. These beds can be effortlessly wheeled to a different place thanks to their castors. This enables access for care from all sides as well as temporary spatial separation of the partners if this is required for medical reasons.

### Cosy and versatile

The Regia partner is available in many attractive wood decors and offers upholstery options with fabric or faux leather. It lends itself to every interior. The Easy Switch system enables split safety sides to be fitted without tools at any time. Two handsets are available for adjusting the two mattress bases.

### Features

Mattress base: 2 x 90 x 200 cm, in 4 sections, with easy-clean metal slats; mattress base height adjustment from 26 – 80 cm, all mattress base elements electrically adjustable, reverse-Trendelenburg position up to 15° with automatic auto contour setting; mattress base with automatic backrest retraction; bed frames move synchronously during height adjustment and tilting for ideal pinch protection; safe working load 2 x 225 kg


\*optional accessories



▲ Nobla (wood decor or upholstered, longer delivery time)



▲ Nobla at head end, Nobla mala at foot end (wood decor or upholstered)

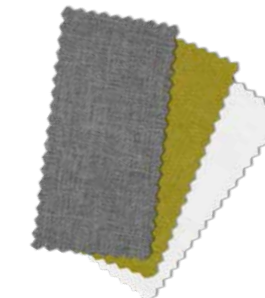


▲ Fabric: Light green Tiba 6210; faux leather: White 1-1142-197, Anthracite 1-1142-395



▲ Natural Beech



▲ Havana Cherry

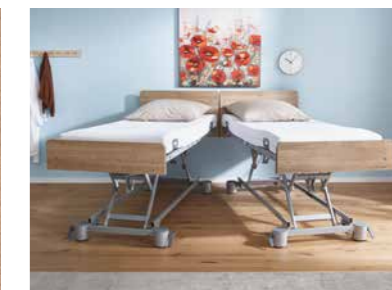

▲ Separable bed halves



▲ One handset each per mattress base



▲ The comfort mattress base comprises 50 free-floating spring elements per bed half.\*



▲ The pivoting standing aid makes it easier to get up out of bed in an ergonomic way.



*Regia gardo*



---

## *Regia gardo* The bed for optimum all-round protection

---

### Contemporary solutions for particular requirements

Diseases such as Huntington's disease, whose symptoms include uncontrolled movements, expose persons to a higher risk of impact injuries. Their beds are often individually padded, for example with foam pieces cut to size. The Regia gardo offers a contemporary solution that combines safety, comfort, and optimal hygiene.

### Tailor-made upholstery

The Regia gardo combines all the strengths of a modern electrically-adjustable care bed with custom-made padding on the inside. The headboard, footboard and full-length safety sides are equipped with padded covers that are easy to remove for thorough cleaning. The bed has an attractive overall appearance and is pleasing to the touch, which makes it a place of dignified peace and relaxation.

### Safe working load of 225 kg

The Regia gardo also offers enhanced comfort for heavier residents.

### Features

Mattress base: 90 x 200 cm, in 4 sections, with sturdy, easy-clean metal slats; mattress base height adjustment from 26 – 80 cm, all mattress base elements electrically adjustable; reverse-Trendelenburg position up to 15° with automatic auto contour; mattress base with automatic backrest retraction; removable padding on the interior of the bed; full-length safety side; Nobla headboard and footboard



▲ The padded covers are fitted closely to the headboard and footboard, thus preventing the creation of unhygienic gaps.



▲ Segufix locks keep the ends of lateral covers positioned tightly against the bed frame.



▲ The pieces of padding can be easily removed for thorough cleaning.



▲ This seamless all-round protection keeps residents safe from impact injuries.



▲ The optional wall deflection roller prevents collision damage while moving the bed and adjusting its height.



▲ Convenient Velcro straps ensure that no gaps appear between the lateral covers and the bed frame.



*Inovia*

# Inovia

## The care bed packed with innovative ideas

### Low-height bed for fall prevention

The height adjustment range featured by the Inovia extends from a very low 22 cm to 77 cm. This is low enough to allow efficient fall prevention as well to provide a convenient height for the resident to get in and out of bed. It also provides an ergonomic waist level height for caregivers to work at.

### Telescopic safety sides

The optional telescopic safety sides (TSG) are easy to raise in 2 stages. If they are only raised at the head end, this leaves ample room for getting in and out of bed.

### Safety side protective height of 43 cm

Mattresses and anti-decubitus systems up to a maximum thickness of 21 cm can be used in the care bed without the need for adjustable safety side height extensions.

### 24-volt system included as standard

The 24-volt drive system dispenses with the requirement for leakage current measurements by the operator of the bed during the expected service life of the bed and beyond. The length of this lifetime is given in the operating instructions. Residents are provided with optimal protection as all live 230 volt components are kept outside the care bed.

### Safe working load of 225 kg

The care bed is also suitable for heavier residents.

### Features

The Inovia is also available with a mattress base width of 100 cm.

Mattress base 90 x 200 cm or 100 x 200 cm, in 4 sections, backrest and thigh rest and mattress base height electrically adjustable; castors lockable in pairs in every mattress base position; mattress base with automatic backrest retraction; wooden surround in beech decor/ solid beech, optionally also in Real Cherry, Havana Cherry or Royal Maple

\*optional accessories



▲ Solid wood handle bar



▲ Castors can be locked in pairs in every mattress base position



▲ Backrest retraction



▲ Handset with selective locking function



▲ Comfort mattress base\*



▲ Tray\*



▲ Reading lamp\*



▲ Telescopic safety side\*



▲ Height adjustment range: 22 cm to 77 cm



*Westfalia-Klassik*

# Westfalia-Klassik A feeling of home whatever your style

## A bed for the discerning

The Westfalia Klassik is designed for discerning users who are not prepared to compromise on comfort, individuality and sophisticated technical features in their daily care routine. This care bed looks good in every situation: If the full-length safety side is lowered, its bars rest seamlessly on top of each other. If it is raised, attractive wooden panels cover the bed frame along both sides of the bed.

## Height adjustment range from 35 to 80 cm

With its adjustable height range of 35 to 80 cm, the Westfalia Klassik contributes to fall prevention, can be easily adjusted to the ideal height for getting in and out of bed and allows caregivers to work at a back-friendly height.

## Features

Mattress base 90 x 200 cm, in 4 sections, backrest, thigh rest and mattress base height electrically adjustable; automatic backrest retraction; attractive wooden surround in Natural Beech with integrated full-length safety sides; framework colour Topaz; metal mattress base; central castor locking; safe working load 200 kg; 24-volt drive system



▲ Central locking castors



▲ Electronic switch mode power supply unit integrated immediately downstream of the power plug



▲ Mattress base height adjustable from 35 to 80 cm, adjustment range = 45 cm



▲ Wooden surrounds: Standard, Comodo with round posts\*, Primero with handle bar\*



▲ Tray\*



▲ Wall deflection roller\*



▲ Reading lamp\*



▲ Bed extension with metal support bars, safety sides and side panels\*



▲ Foam leather cover for safety sides\*

\*optional accessories

Mattress  
base width:  
90 cm  
100 cm



*Dali with wooden surrounds*

# Dali with wooden surrounds

The practical care bed that brings a touch of cosiness

## Quick and dependable

The Dali uses the advantages of digital technology to make life easier for residents and their families. In addition to the cabled handset, the wireless model is convenient to use and offers great freedom of movement. Its adjustment functions can be conveniently locked and enabled on both handsets. Dali is also available on request with an environmental control, such as speech or eye control.

## Attractive wooden surrounds in “Arminia style”

Headboards and footboards on the Dali can be optionally covered with wooden surrounds. This attractive feature on the Dali takes up an idea from the previous model Arminia.

## Safety sides with Easy Click system

With the Easy Click system, the full-length safety side can be fitted and removed with a single click, without the need for tools. This makes it easier to set up the bed and offers flexibility in the care situation.

## Features

Mattress base 90 x 200 cm or 100 x 200 cm, in 4 sections, height adjustable from 40 – 80 cm, backrest, thigh rest and mattress base height electrically adjustable; reverse-Trendelenburg position adjustable by tilting the mattress base; handset with selective locking function; safe working load 185 kg

## Stylish alternative: Dali with Softcovers

If you prefer a bed with a fabric cover rather than a wooden surround, our stylish Softcovers offer an interesting alternative. These can be used to cover the headboard and footboard as well as the safety sides of the Dali. The machine-washable fabrics in attractive colours transform the bed into an individual designer piece.

\*optional accessories



▲ Reverse-Trendelenburg position



▲ Convenient standing aid instead of full-length safety side\*



▲ Cabled and wireless handset



▲ Adaptable safety sides\*



▲ Tray\*



▲ Bed extension including safety side bars\* (also for retrofitting)



▲ Reading lamp\*



▲ Foam leather cover\*



▲ Natural Beech



*Dali economic*



# Dali economic

## The new generation in home care

The Dali economic offers individual options in everyday care.

### basic variant

- The budget-oriented bed without a safety side and thigh-rest adjustment for independent residents with low protection needs
- Offers health-promoting reverse-Trendelenburg position

### pro care variant


- The comfortable variant offering thigh rest adjustment, a wooden safety side and a patient lifting pole

### Features


Mattress base 90 x 200 cm, in 4 sections, mattress base height adjustment of 40 – 80 cm, backrest and mattress base height adjustable with electric motors, thigh rest adjustment and reverse-Trendelenburg position optional; handset with selective locking function; safe working load 185 kg




▲ The wooden slat mattress base is included in both standard stock variants, while the full-length wooden safety side is included in the pro care variant.

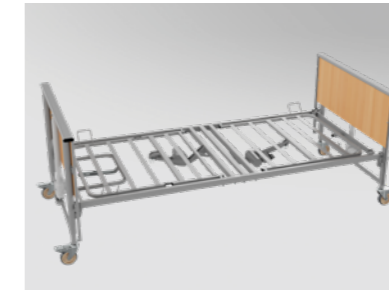

▲ Metal mattress base\*

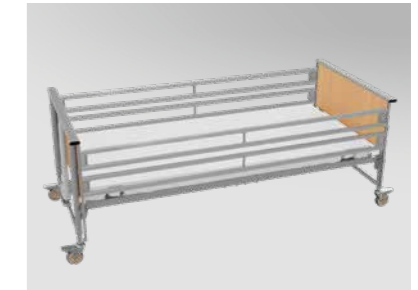

▲ Metal safety sides\*

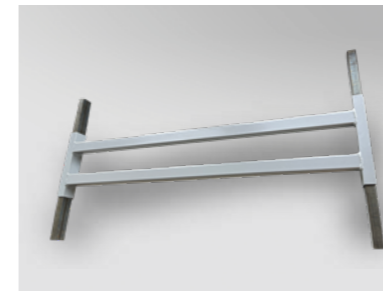

▲ Bed extension\*



▲ Two piece lifting pole with triangular grab handle



▲ Wooden surround\*




▲ Ratchet mechanism with Fowler position\*

Our easy to use handsets are adapted to the functions offered by the variants. Additional functions and suitable handsets can be individually configured.



◀ HL 73 “Standard” is included in the basic version

HL 73 “pro care” is included in the pro care variant ▶

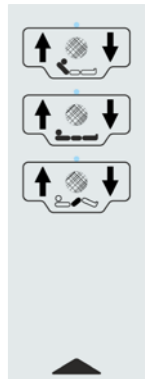

\*optional accessories



*Lippe IV bed insert*

# Lippe IV bed insert

For fitting to an existing bed

## Comfort in your familiar bed

Many people would like to enjoy the benefits of an electrically adjustable bed but would prefer to retain their familiar bedroom furniture. Our Lippe range bed inserts can simply be placed inside an existing bed frame. This way, they provide ample support but remain as invisible as a normal slatted frame.

## 24-volt system included as standard

The 24-volt drive system dispenses with the requirement for leakage current measurements by the operator of the bed during the expected service life of the bed and beyond. The length of this lifetime is given in the operating

instructions. Residents are provided with optimum protection as all live 230 volt components are kept outside the bed insert frame.

## Features

Mattress base in 4 sections, backrest, thigh rest and mattress base height adjustable electrically; lower leg rest manually lowerable with a ratchet, mattress base adjustment from 29 – 74 cm; chassis adjustable individually to slatted bed frame widths of 100 and 90 cm; 5th mattress retainer bar included as standard for holding the mattress in place at the foot end: safe working load 200 kg; handset with selective locking function

## External dimensions

Nominal mattress base size	Actual mattress base size
80 x 190 cm	75 x 185 cm
80 x 200 cm	75 x 195 cm
90 x 190 cm	85 x 185 cm
90 x 200 cm	85 x 195 cm
100 x 190 cm	95 x 185 cm
100 x 200 cm	95 x 195 cm

\*optional accessories



▲ Mattress base extension of 10 cm for inserting at head end or foot end.



▲ Handset with selective locking function



▲ Foot with floor-level adjustment screw and fixing strap



▲ Power socket for connecting the electronic switch mode power supply unit



▲ Patient lifting pole\*



▲ Safety sides\* (also for retrofitting)



▲ Standing aid\*



▲ Adjustable feet for raising the lowest position (+ 7 cm)



▲ Widening rail, 5 cm\*



*Lippe 120/140 bed insert*

# Lippe 120/140 bed insert

The height-adjustable bed insert with a comfortable width

## Quality and longevity

The Lippe 120 and Lippe 140 bed inserts are extremely sturdily built and are ideal for meeting the challenges of day-to-day care. Their high quality ensures a long useful life with low maintenance and decisively minimises follow-on costs. Mattress base widths of 120 or 140 cm meet all discerning expectations when it comes to maximising comfort.

## 24-volt system included as standard

The 24-volt drive system dispenses with the requirement for leakage current measurements by the operator of the bed during the expected service life of the bed and beyond. The length of this lifetime is given in the operating instructions. Residents are provided with optimum protection as all live 230 volt components are kept outside the bed insert frame.

## Safe working load of 250 kg

The bed insert is also suitable for heavier residents.

## Features

Mattress base in 4 sections, backrest, thigh rest and mattress base height adjustable electrically; lower leg rest manually lowerable with a ratchet, mattress base adjustment from 29 – 74 cm; chassis with feet including floor-level adjustment screws and 4 straps for fixing within the existing bed; 5th mattress retainer bar included as standard for holding the mattress in place at the foot end; handset with selective locking function; optional 20 cm bed extension, also for retrofitting

### External dimensions

	Lippe 120	Lippe 140
Nominal mattress base size	120 x 200 cm	140 x 200 cm
Actual mattress base size	115 x 195 cm	135 x 195 cm

\*optional accessories



▲ Foot with floor-level adjustment screw and fixing strap



▲ Ball bearing height adjustment



▲ Electronic switch mode power supply unit integrated immediately downstream of the mains plug



▲ Power socket for connecting the electronic switch mode power supply unit



▲ Large height adjustment range of 29-74 cm



▲ Standing aid\*



▲ Adjustable feet for raising the lowest position\* (+7 cm)



▲ Bed extension\* (also for retrofitting)



▲ Patient lifting pole\*



*Relax*

# Relax

## The combinable bed frame

### Stylish bed surround for a contemporary dream bed

Lippe bed inserts fit in practically every familiar bed – but sometimes giving the bed a new look as well enhances the feeling of well-being. Our stylish contemporary Relax is the ideal bed surround for combining with a Lippe or an electric slatted frame. It is available in a choice of 6 sizes: 90 x 200 cm, 100 x 200 cm, 120 x 200 cm, 120 x 220 cm, 140 x 200 cm and 140 x 220 cm. Two Relax beds can also be joined together easily with bolts to make a double bed.

### Decors to suit individual tastes

The Relax is available in four different wood decors, which offer a choice for every furnishing style. Enjoy the cosy atmosphere exuded by the pale gold of the Royal Maple, the subtle Beech, the smoky Havana Cherry or the rich Bella Noce Choco. If you use your Relax in combination with a Lippe bed insert, you can conceal its feet behind a stylish grey bed surround.

### Easy to get in and out of bed

The top side edge of the bed is 45 cm above the floor – that is an ideal height for getting in and out of bed for many people in need of care. The use of height-adjustable strips for supporting the slatted frame makes it possible to adjust this height even more closely to individual needs. With a Lippe bed insert, the mattress base can be raised to 74 cm, which means that caregivers can work at a back-friendly ergonomic height.

### Accessories

When combined with a Lippe IV or Lippe 120/140 bed insert, it goes without saying that the Relax can take advantage of the accessories available for these bed inserts – such as the patient lifting pole or the standing aid. In addition, customers can benefit from the safe and cost-saving 24-volt drive system for the Lippe models.

\*optional accessories



▲ Slatted frame electrically adjustable.




▲ Floor-level bed sides.




▲ Safety sides\* (also for retrofitting).




▲ Standing aid in combination with a bed insert from the Lippe range\*.

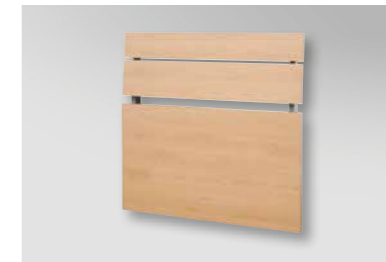

▲ Wooden frame in Natural Beech; panels and bed sides in Natural Beech decor.

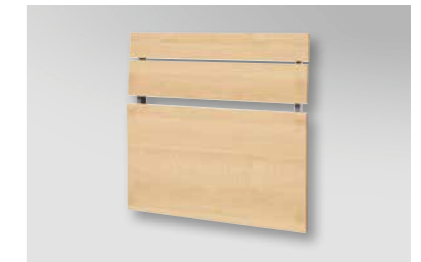

▲ Wooden frame in Royal Maple; panels and bed sides in Royal Maple decor.

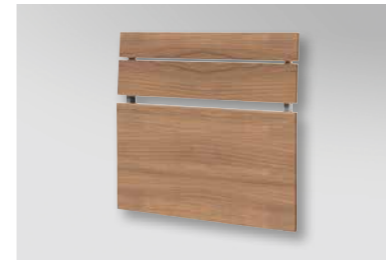

▲ Wooden frame in Havana Cherry decor; panels and bed sides in Havana Cherry decor.



▲ Wooden frame in Bella Noce Choco decor; panels and bed sides in Bella Noce Choco decor.



*Gigant*



# Gigant The sturdy heavy-duty care bed

## Quality and longevity

The Gigant care bed is extremely sturdily constructed and stands up to all the rigours of daily care. Their high quality ensures a long useful life with low maintenance and decisively minimises follow-on costs. People in need of care will feel right at home in this bed.

## Safe working load of 350 kg

The mattress base can be loaded with a weight of up to 350 kg, including the mattress and accessories. Residents weighing up to 280 kg will appreciate the ideal and generous 120 cm-wide mattress base. An optional bed extension to 220 cm is available to choose from.

## Mobility and comfort

Mobilisation of residents is becoming increasingly important in the daily care routine. The mobilisation support for the optional split safety side for the Gigant is an ideal aid for standing up due to its ergonomic knob. Raised only at the head end, the split safety side also creates a good balance between safety and freedom for the resident.

## Features

Mattress base 120 x 200 cm, one-piece, not dismantlable, mattress base height adjustable from 40 – 80 cm, backrest, thigh rest and mattress base height are adjustable with electric motors; reverse-Trendelenburg position adjustable by pivoting the mattress base using the hand-set; emergency backrest release; metal parts sealed with polyester powder coating (colour Argentum), wooden surround in Beech decor/solid Beech

\*optional accessories



▲ Mobilisation aid holder (left) and locking box (right)



▲ Emergency backrest release



▲ Extremely sturdy safety side construction



▲ KST panels on the chassis discreetly conceal the technical components



▲ Clamp-on extra-high safety side bar\*



▲ Clamp-on extra-high safety side bar\*



▲ Reading lamp\*




▲ Bed extension, including longer safety sides and side panels\*



▲ Foam leather cover\*

# Mattresses

More information and all available dimensions at:  
[burmeier.com/mattresses](http://burmeier.com/mattresses)

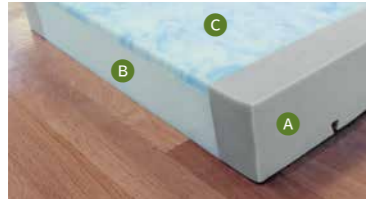

	Com- pression hardness	Density
A	5.0 kPa	50 kg/m <sup>3</sup>
B	4.0 kPa	
C	2.5 kPa	40 kg/m <sup>3</sup>

## Sky-fit – the premium model

- High-quality dual-layer mattress with a 2 cm-Skytherm® top layer and a 12 cm cold foam base layer
- Incisions in the base layer ensure it optimally adapts to the mattress base
- Reinforced edge zone allows users to sit comfortably on the side of the bed
- Removable reversible cover, one side polyurethane (PU), the other side polyester and viscose, which can be washed in a standard washing machine



	Com- pression hardness	Density
A	5.0 kPa	-
B	4.0 kPa	40 kg/m <sup>3</sup>

## Comfort-fit – the mid-range mattress

- The comfortable mattress with a core made of resilient polyurethane cold foam (14 cm)
- Incisions in the base layer ensure it optimally adapts to the mattress base
- Edge zone reinforcement protects against falling out of bed
- Removable reversible cover, one side polyurethane (PU), the other side polyester and viscose, which can be washed in a standard washing machine



	Com- pression hardness	Density
A	4.5 kPa	40 kg/m <sup>3</sup>

## Basic-fit – the entry-level model

- Single-layer mattress made of cosy polyurethane (12 cm)
- White removable PU cover, can be washed in a standard washing machine



	Com- pression hardness	Density
A	5.0 kPa	50 kg/m <sup>3</sup>
B	9.0 kPa	60 kg/m <sup>3</sup>
C	2.5 kPa	40 kg/m <sup>3</sup>

## Cloud-fit – the versatile heavy-duty mattress

- Particularly suitable for very heavy residents with a body weight of 140 to 240 kg
- 3-layer mattress with a top layer made of comfortable Sky-therm® foam, and two resilient cold-foam layers



*Bedside  
cabinets*

## Bedside cabinets and side tables



### Cherusker bedside cabinet

With height-adjustable tilting overbed table, plastic-coated body, beech decor

#### Technical data:

Height: 87 cm  
Width: 52 cm  
Overall depth: 50 cm  
Overall width (overbed table horizontal): 114 cm  
Overbed table height adjusts from 82 to 112 cm  
Double castors Ø 50 mm, 2 locking castors



### Arminius bedside cabinet

With height-adjustable tilting overbed table, plastic-coated body, beech decor

#### Technical data:

Height: 89 cm  
Width: 54.5 cm  
Overall depth: 45 cm  
Overall width (overbed table horizontal): 125 cm  
Overbed table height adjustment from 78 to 112 cm  
Double castors Ø 50 mm, 2 locking castors



### Hermann bedside cabinet with overbed table

With fold-down overbed table, plastic-coated body, beech decor

#### Technical data:

Height: 89 cm  
Width: 52.5 cm  
Overall depth: 45.5 cm  
Overall width (overbed table horizontal): 114 cm  
Double castors Ø 50 mm, 2 locking castors



### Varus bedside cabinet

Plastic-coated body, no overbed table, Available in four different decors: Natural Beech, Royal Maple, Havana Cherry, Bella Noce Choco

#### Technical data:

Height: 57 cm  
Width: 54 cm  
Overall depth: 45 cm  
Double castors Ø 50 mm, 2 drawers



▲ Duo model



▲ Solo model

### Overbed tables and side tables Duo and Solo

With height-adjustable, tilting overbed table (Duo: additional fixed side tray), Beech decor

#### Technical data:

Height: 81 cm  
Width: 60 cm (Solo), 81 cm (Duo)  
Overall depth: 42 cm  
Overbed table height adjusts from 80 to 113 cm  
Castors Ø 50 mm, 2 locking castors



### Pleto server

With height-adjustable tilting overbed table and fixed side tray, available in 4 decors: Natural Beech, Royal Maple, Lindberg Oak and Havana Cherry.

#### Technical data:

Width: 90.2 cm  
Height: 69 – 104 cm  
Overall depth: 42 cm  
Tilting overbed table 30° – 60°  
Castors Ø 50 mm, 2 locking castors

For your notes



# STIEGELMEYERGROUP

## Solutions for people who care

The Stiegelmeier Group is Germany's leading manufacturer of medical beds for hospitals, nursing homes and homecare. With modern bedside cabinets, furniture and digital assistance systems, it is also an innovative provider of comprehensive solutions. The family-run company successfully develops offers for international markets in more than 100 countries.

### Our strengths:

- Product development at our Herford headquarters in Germany – production in four modern factories in Germany and Poland
- Factory in China for the Chinese market
- Subsidiaries and sales partners in many countries
- Top-class service with reliable spare parts supply
- Training and process optimisation for our customers
- Sustainable action with strong social commitment



1 000 EMPLOYEES WORLDWIDE	11 LOCATIONS in Germany, Poland, the Netherlands, Belgium, France, Finland, South Africa and China	 Tailor-made <b>OFFERS</b> for all <b>MARKETS</b>
MORE THAN 100 000 beds produced every year		International <b>CERTIFI- CATION</b>
UP TO 15-YEAR spare part availability	<b>DISTRIBUTION</b> in more than 100 countries	 <b>GERMAN ENGINEERING</b>

### Our expertise

- Hospital
- Nursing home
- Homecare
- Lifestyle beds
- Furniture systems
- Service

### Our values

- Experienced
- Human
- Reliable
- Solution-focused
- Open-minded
- Confident

### Our brands





# STIEGELMEYERGROUP

Solutions for people who care

▼ Your specialist dealer



## Hauptsitz / Headquarters

*Deutschland / Germany*  
**Burmeier GmbH & Co. KG**  
Industriestraße 53  
32120 Hiddenhausen  
Phone +49 (0) 5223 9769 - 0  
Fax +49 (0) 5223 9769 - 090  
info@burmeier.com  
www.burmeier.com



## Hauptsitz / Headquarters

*Deutschland / Germany*  
**Stiegelmeier GmbH & Co. KG**  
Ackerstraße 42  
32051 Herford  
Phone +49 (0) 5221 185 - 0  
Fax +49 (0) 5221 185 - 100  
info@stiegelmeier.com  
www.stiegelmeier.com

Delivery via medical equipment retailers only

### *France*

**Stiegelmeier SAS**  
16 Avenue de l'Europe  
Bâtiment SXB1  
67300 Schiltigheim  
Phone +33 (0) 4 7238 - 5336  
contact@stiegelmeier.fr  
www.stiegelmeier.fr

### *Nederland / Netherlands*

**Stiegelmeier BV**  
Peppelenbos 7  
6662 WB Elst  
Phone +31 (0) 345 581 - 881  
info@stiegelmeier.nl  
www.stiegelmeier.nl

### *Polska / Poland*

**Stiegelmeier Sp. z o.o.**  
86-212 Stolno  
Grubno 63  
Phone +48 56 6771 - 400  
Fax +48 56 6771 - 410  
info.stolno@stiegelmeier.pl  
www.stiegelmeier.pl

### *South Africa*


**Stiegelmeier Africa (PTY) Ltd.**  
The Estuaries - Century City  
Cape Town, 7441  
Phone +27 (0) 21 555 - 1415  
Fax +27 (0) 21 555 - 0639  
info.southafrica@stiegelmeier.com  
www.stiegelmeier.co.za

### *Suomi / Finland*

**Stiegelmeier Oy**  
Ahjokatu 4 A  
15800 Lahti  
Phone +358 (0) 3 7536 - 320  
stiegelmeier@stiegelmeier.fi  
www.stiegelmeier.fi

### *China*

**Stiegelmeier China Co. Ltd.**  
Poly Center, No. 8 South Wenhua Road  
Dongping New Town, Shunde District  
528300 Foshan  
Phone +86 757 29878710 - 0  
info@stiegelmeier.cn  
www.stiegelmeier.cn

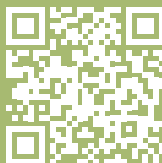

stiegelmeier-group.com

Nicht alle Produkte und Konfigurationen sind in jedem Land verfügbar. Not all products and configurations are available in every country.  
Technische Änderungen vorbehalten. Leichte Farbabweichungen der Originale sind drucktechnisch bedingt. Stand 11/2024. Version 14, Artikel-Nr. 601725.  
Subject to technical changes without notice. The shown product colours may differ slightly from the actual product. Last updated 11/2024. Version 14, Item no. 601725.