


SAFE LIVING AT HOME

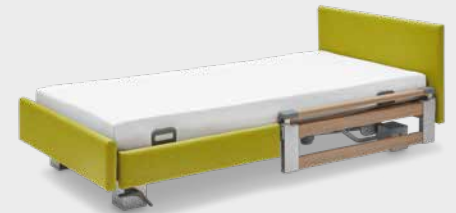

Retailer catalogue

Beds and accessories for caring at home

Ready for
the MDR from
May 2020

International Edition

STIEGELMEYERGROUP



Contents

General information


5	Burmeier, the company	10	Econ drive system
6	The Stiegmeyer-Group	12	Easy Switch
8	Safety standards for care beds	14	The Burmeier apps
9	The 24-volt drive system	16	Technical data

Care beds and bed inserts

24	Dali	60	Regia
28	Dali wash	66	Inovia
32	Dali low-econ	72	Westfalia IV
36	Dali low-entry	76	Westfalia-Klassik
42	Lippe IV	82	Gigant
46	Lippe IV washable	86	Allura II
50	Lippe 120/140	92	Vertica homecare
54	Relax		

Bedside cabinets and furniture

98	Bedside cabinets and side tables	110	Wardrobes and cabinets
104	Ravello	113	Decor, tables
108	Mattresses	114	Chairs

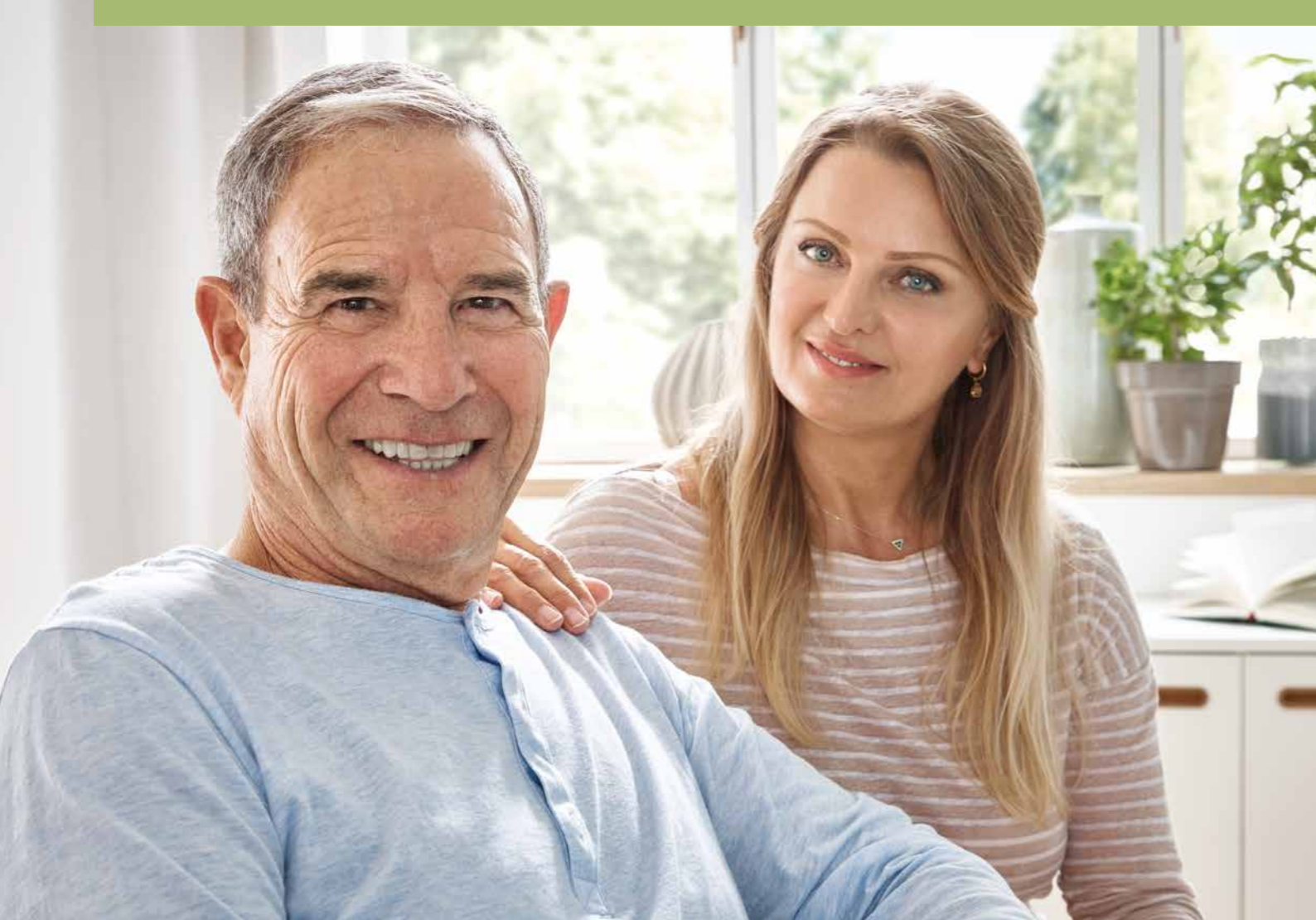

Safe living at home

Burmeier is your specialist for homecare beds and furniture – and has been the market leader in Germany for many years. Our aim is to make life at home as easy and pleasant as possible for those in need of care and for their caring relatives and to promote independence and mobility. Our products are technically innovative and durable, and enjoy an excellent reputation worldwide thanks to their high quality and safety standards.

In its current form, Burmeier has been another successful subsidiary of the Stieglmeyer-Group since 1995. Together, the companies in the group can look back on more than 120 years' experience in manufacturing and selling hospital and care beds and accessories. The combined innovative force of the corporate group focuses on developing innovative technologies to shape the future in all care sectors. Examples include digital assistance systems and a range of safety side concepts to give greater safety and freedom to those in need of care and also sustain mobility.

The comfortable ambience, high functionality, and fast availability of Burmeier products are decisive arguments for medical product retailers. Factors such as virtually constant availability of all articles, large stocks of products, and a reliable logistics concept ensure that products can be delivered to anywhere in Germany at only 24 to 48 hours' notice. Burmeier uses identical parts in a variety of models – creating sustainable added value for retailers when managing and organising procurement.

The right bed for every case: this is the aim of the company group, in a nutshell. This means that an affordable entry model bed can be found in this catalogue alongside heavy-duty beds and low-height beds. Discover in this catalogue how Burmeier products can support all those caring for others at home and can make life easier and more comfortable.



Custodians of values Shapers of the future

Custodians of values

Our company was founded in Herford in 1900 and is still a family-run enterprise today. Respect and a sense of responsibility towards our customers, employees and partners form the basis for our day-to-day work. Based on these values and our many years of experience, we offer our customers all over the world high-quality products, solutions and services to meet the needs of the care and healthcare sectors.

Shapers of the future

In order to cope with rapid changes in the care and healthcare sectors, we develop all our products ourselves. Creating innovations and visions of what care will be like in future are the goals we pursue. We maintain regular contact with our customers in order to support them in their daily tasks. Implementing tomorrow's trends today is a successful tradition of ours.

Our expertise

Hospital

Nursing home

Homecare

Lifestyle beds

Furniture systems



The Stiegmeyer-Group is your international partner for hospitals and care facilities and, with the Burmeier brand, also for caring in the home. We support our customers with practice-oriented products, solutions and services that we are constantly developing further.

- Founded in 1900 and still family managed
- Over 1,100 employees worldwide
- All core products manufactured in proprietary factories
- Product development and production in Germany
- Service, training and process optimisation at the customer's premises
- International certification and export to more than 60 countries
- Distribution companies located in Belgium, Finland, France, Netherlands, Poland and South Africa

Our brands



Our key focuses and innovations



Safety standards

All the beds are designed for maximum safety as standard.

This includes:

- Primary fuse, electronic overload cut-off device and thermal fuse against overheating
- Motors run on safety extra-low voltage
- In-company quality standards exceed industry standard requirements
- Anti-kink strain relief on bed frame and on motor housing
- At least IPX4 splash-proof
- Optional IPX6 splash protection for automated reprocessing

The following standards and directives are applied:

- Up to May 2020: Medical Devices Directive 93/42/EEC
- From May 2020 onwards: Medical Device Regulation (MDR) (2017/745)
- DIN EN 62353: Medical electrical equipment – Recurrent test and test after repair of medical electrical equipment
- DIN EN 12182: Technical aids for disabled persons
- DIN 32977: Adapted designing for disabled persons
- DIN EN 60601-1: Medical electrical equipment – General requirements for basic safety and essential performance
- DIN EN 60601-1-2: Medical electrical equipment – Collateral standard: Electromagnetic compatibility
- DIN EN 60601-1-6: Medical electrical equipment – Collateral standard: Suitability for intended use
- DIN EN 60601-2-52: Medical electrical equipment – Particular requirements for basic safety and essential performance of medical beds
- DIN EN ISO 14971: Medical devices – Application of risk management to medical devices



24-volt drive system


Burmeier's innovative 24-volt drive system provides an ideal safety standard for caring for people at home. All live 230-volt components are banished from the care bed. Immediately downstream of the mains plug is a housing containing the **SwitchMode** electronic voltage transformer. Residents therefore enjoy a high level of safety.

The bed has no current-carrying components except when the handset is in use. Power consumption drops by up to 93% as a result. Even more costs and effort are saved by the fact that **leakage current measurements are no longer necessary**, as certified by **TÜV Rheinland**. In principle, the required leakage current measurement must be carried out every 2 years for care beds in commercial use. The use of the 24-volt system prolongs this interval so that measurements are not necessary within the expected lifetime of the bed. The applicable lifetime is indicated in the operating instructions.

Even when taking the bed into service, leakage current measurements are unnecessary as these have already been carried out in the factory. It is also not necessary for **an electrically trained person** to be involved in assembling the bed – this saves the bed operator additional costs.

SAFE LIVING AT HOME


Energy saving

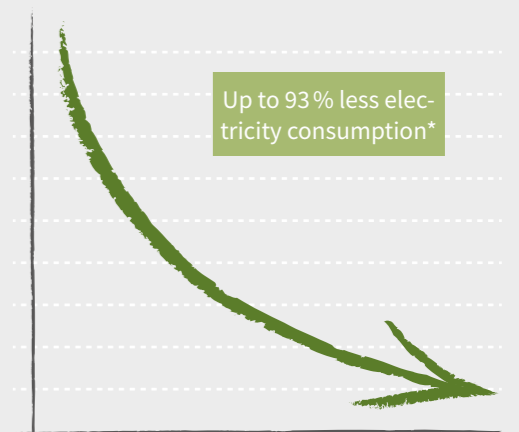

* Compared to 230-volt systems, assuming approx. 5 adjustments per day

SAFE LIVING AT HOME

Strong advantages

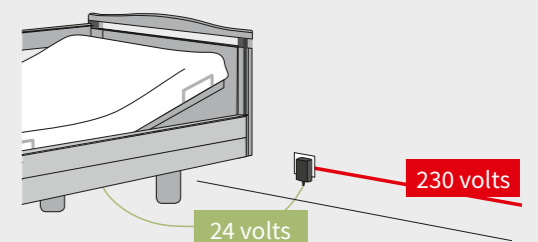

- Safety extra-low voltage starting from the transformer at the power socket
- No 230-volt components on the bed
- No leakage current measurements are necessary after installation or beyond during the expected service life of the bed
- A test certificate can be downloaded for every bed.



For the models:
Dali econ
Dali low-econ
Dali low-entry-econ



Econ drive system



▲ The budget-oriented low-height bed Dali low-econ.



▲ The low-height bed Dali low-entry-econ for perfect fall prevention.



▲ The Econ drive system is available with many different international power plugs.

Most bed models in the Dali range are equipped with an electricity-saving 24-volt drive system as standard. In some countries, however, power plug standards apply that are not compatible with the 24-volt system and its SwitchMode transformer unit plug unless an accessory is used. In this case, our Econ models with an **alternative drive** are an attractive and cost-effective choice.

Both the Dali and the Dali low-entry are available in the versions **Dali econ** and **Dali low-entry-econ**. A third model, the **Dali low-econ**, is also available, which does not exist in a 24-volt version.

All Econ beds have a safe, reliable and economical **230-volt drive**. With a selection of suitable power plugs to choose from, these care beds can be used internationally, wherever they are needed. The **handsets** of the Econ models are connected to the bed with a cable, as was previously common practice. The handset can be used to selectively lock and unlock all adjustment functions. The Dali-Lock app cannot be used with the Econ beds.



▲ The corded handset of the Econ models can be fixed anywhere on the safety sides.



Easy Switch for our Regia bed

The patented system for adaptable care beds with standard-compliant safety


Choose the Regia with the telescopic split safety side or without a safety side and benefit from our flexible Easy Switch system. The Easy Switch principle is as easy as it is impressive: The safety side elements, the headboard and footboard as well as the side panels can be attached and removed **without tools**. This is possible thanks to easy-to-use locks that snap into place securely and cannot be accidentally opened by the resident. Headboards and footboards from the Nobla range and also panels are available as an option with **attractive fabric or faux leather upholstery**.

Thanks to the Easy Switch system, the Regia can easily be fitted any time with safety sides or a new look and **adapted to your requirements**. It becomes a bed that develops with you and supports you in every phase of your life. The removable safety sides are quick to convert to whatever care situation may arise. They also make it considerably easier to assemble the bed.

Modular design for tool-free assembly

Replace headboards in seconds whenever required

- Homelike design variants, with optional handle bars
- Attractive wood decors, fabrics and faux leather to choose from
- Headboard, footboard and side panels freely combinable however you wish

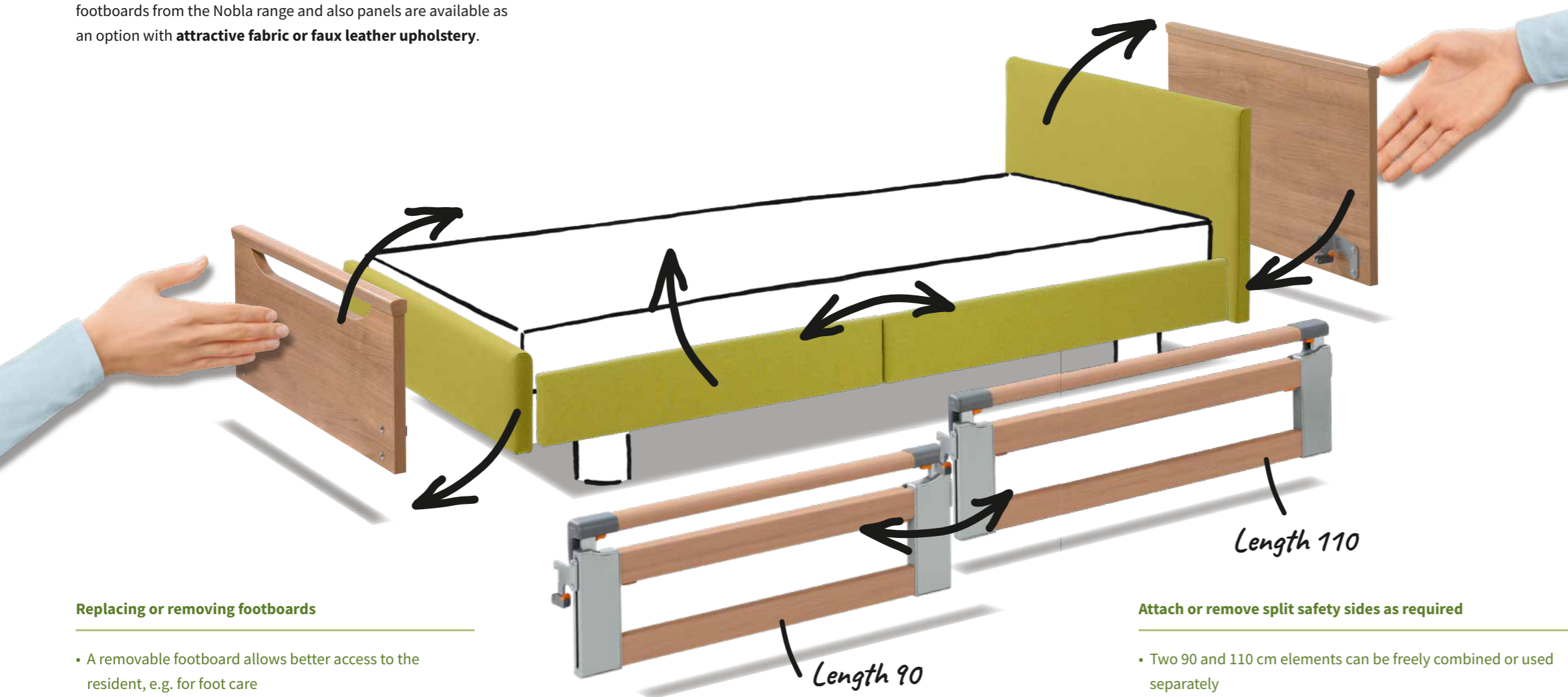

Replacing or removing footboards

- A removable footboard allows better access to the resident, e.g. for foot care
- Choice between a high footboard for a feeling of security or a low footboard for an unobstructed view
- Attractive wood decors, fabrics and faux leather available to choose from

**! Safety sides can be freely combined:
• 110/90 or 90/110**

Attach or remove split safety sides as required

- Two 90 and 110 cm elements can be freely combined or used separately
- Use of only the 110 cm safety side element at the head end provides standard-compliant protection
- Use of only the 90 cm safety side element at the head end provides reassurance and orientation and leaves ample space to get in and out of bed

SAFE LIVING AT HOME

Easy Switch for added flexibility and safety

Beds with Easy Switch meet the different protection needs of residents and can be customised to suit any taste. The flexible safety side provides standard-compliant protection without restricting the mobility of the person in need of care.



▲ All Easy Switch elements are inserted into the bed frame and attached with easy-to-operate locks.



▲ Simply press the fastener upwards.




▲ The lock snaps into place securely – and the element is attached.



▲ To release the fastening again, pull the orange safety catch towards you and press down on the lock.



The future is digital

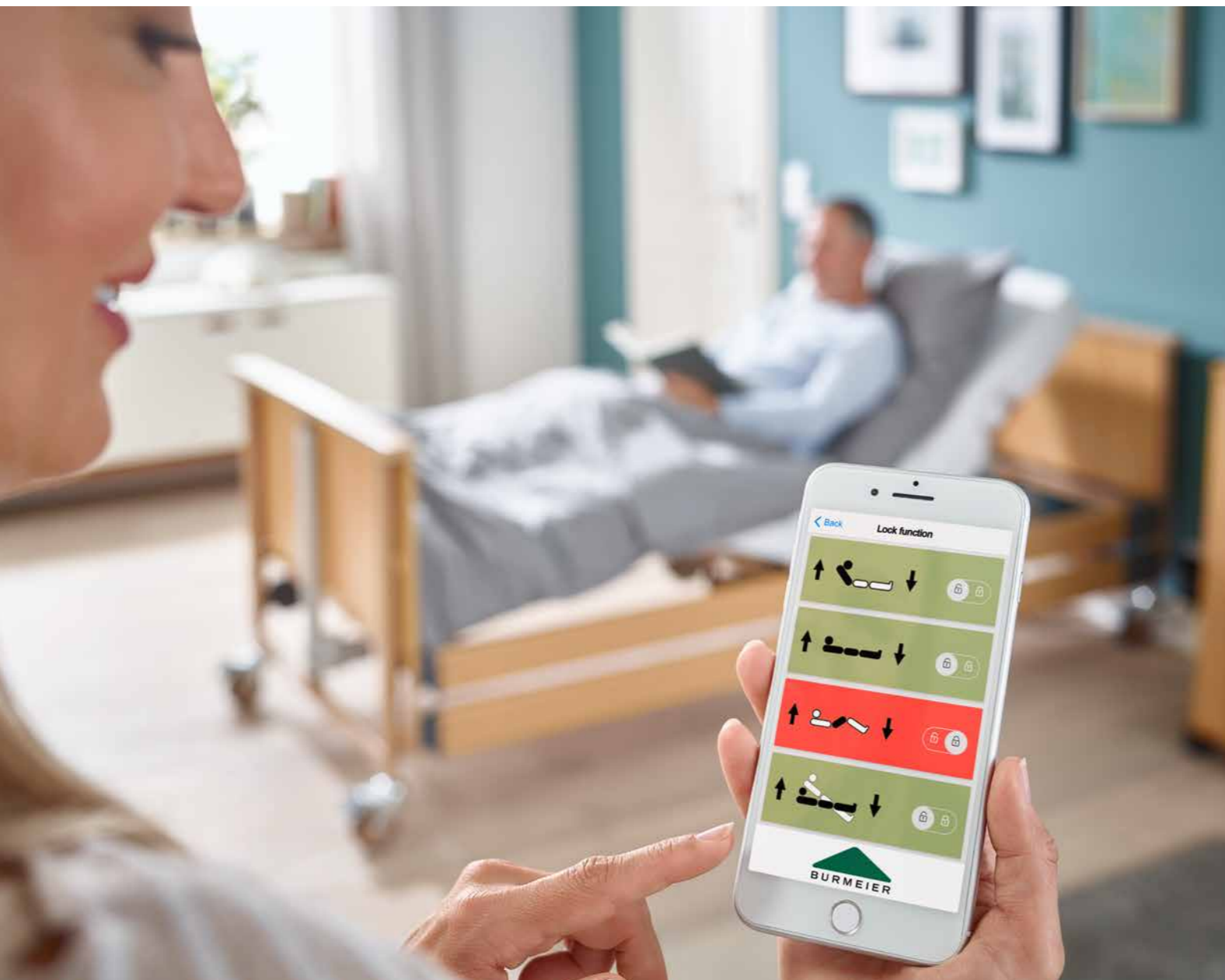

The Burmeier apps

Organise safety digitally

Thanks to the many adjustment options for the mattress base, our Dali beds are comfortable and flexible. Depending on the resident's state of health, however, certain positions may not be advisable. If the carer wants to prevent any use of the reverse-Trendelenburg position or prevent the backrest from being raised, the **Dali-Lock app** is there to assist. This app is available for Android and Apple smartphones and can be downloaded for free from the usual stores.

The Dali-Lock app finds the care bed automatically and connects to it when you scan the code on the back of the handset. Now you can easily block or unblock the individual functions via simple pictograms. The app also keeps you informed of the **battery status** and **activity** of the wireless handset.

The Dali-Lock app is part of Burmeier's **extensive digital network**. Visit our website at www.burmeier.com for news of other apps.

You can find our apps here:



App Store



Google Play




▲ The Dali-Lock app lets you manage the adjustment functions of the bed with ease.



SAFE LIVING AT HOME

The Dali-Lock app

- Download the free Dali-Lock app for your Android or Apple smartphone from the store.
- Start the bed search process. The app will find your Dali.
- Connect your smartphone to the bed by scanning the code on the back of the handset.
- Now you can lock or unlock all backrest or height adjustments of the bed individually.

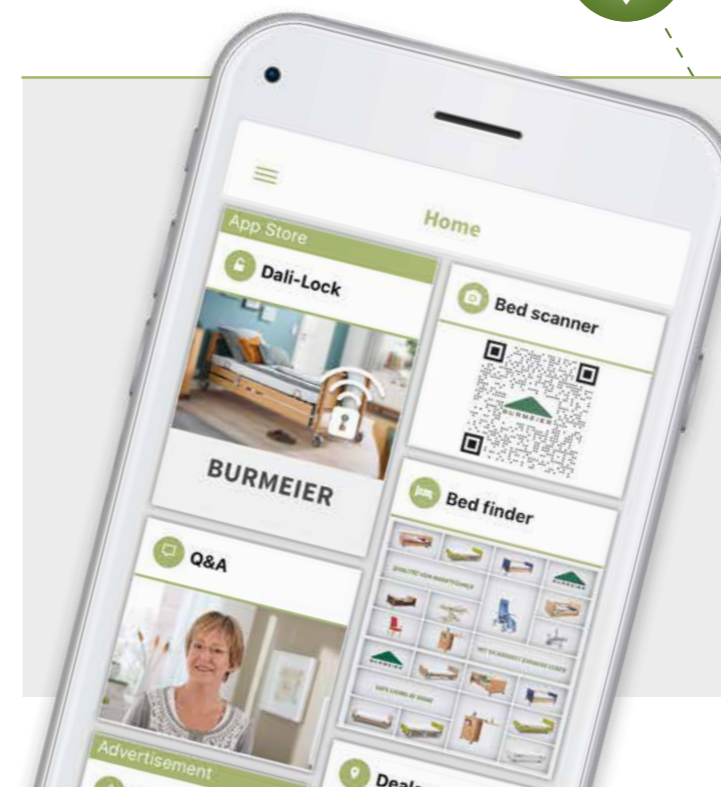

SAFE LIVING AT HOME

The Burmeier app

Download the free Burmeier app from your Apple or Android store.

- Practical bed finder
- Helpful assembly videos
- FORUM customer magazine
- Social media channels
- FAQ

TECHNICAL DATA

Model		Dali	Dali wash	Dali low-econ	Dali low-entry	Lippe IV	Lippe IV washable	Lippe 120/140	Regia
Adjustments									
Height	Electric	✓	✓	✓	✓	✓	✓	✓	✓
Backrest	Electric	✓	✓	✓	✓	✓	✓	✓	✓
Thigh rest	Electric	✓	✓	✓	✓	✓	✓	✓	✓
Locking feature	Selective	✓	✓	✓	✓	✓	✓	✓	✓
Locking feature	Central	–	–	–	–	–	–	–	–
Wireless handset		✓	✓	–	✓	–	–	–	–
24-volt drive system		✓	✓	–	✓	✓	✓	✓	✓
Econ drive system		O	–	✓	O	–	–	–	–
Height adjustment range									
from – to	cm	40-80	40-80	32-72	23-63	29-74	29-74	29-74	26-80
4-section mattress base									
Standard	cm	90 x 200 90 x 180 ¹⁾	90 x 200	90 x 200	90 x 200	90 100 190 200	90 x 190 90 x 200	120 140 x 200	90 x 200
Special mattress base widths		–	–	–	–	80	–	–	100
Metal slats		O	✓	–	✓	–	✓	✓	✓
Wire mesh mattress base		–	–	–	–	–	–	–	–
Sprung wooden slats		✓	–	✓	–	✓	–	–	–
Comfort mattress base		–	–	–	–	–	–	–	O
One-piece		–	–	–	–	–	–	–	✓
Multi-section		✓	✓	✓	✓	✓	✓	✓	–
Adjustable mattress base tilt		✓	✓	✓	✓	–	–	–	✓
Sliding backrest		–	–	–	–	–	–	–	✓
Ground clearance for lifting equipment		✓	✓	✓	✓	–	–	–	✓
Castor locking									
Central		–	–	–	–	–	–	–	✓
Per axle		–	–	–	–	–	–	–	–
Single		✓	✓	✓	✓	–	–	–	–
Castor diameter	mm	100	100	100	75	–	–	–	50
External dimensions and weight									
Length		218	218	225	225	185/195	185/195	195	210
Width		102	102	102	102	85/95	85	115 135	101
Net weight	kg	97	106	97	106	80	80	132 140	135
Safe working load									
	kg	185	185	185	175	200	200	250	225
Wooden surround									
Standard		–	–	–	–	–	–	–	–
Model-related headboard, closed surface		✓	✓		✓	–	–	–	–
Model-related footboard, closed surface		✓	✓	✓	✓	–	–	–	–
Model-related headboard with handle bar		–	–	–	–	–	–	–	–
Model-related footboard with handle bar		–	–	–	–	–	–	–	–
Comodo with round posts		–	–	–	–	–	–	–	–
Primero with handle bar on head and footboard		–	–	–	–	–	–	–	–
Prinzino with handle bar on head and footboard		–	–	–	–	–	–	–	O
Primus closed-surface headboard, footboard with handle bar		–	–	–	–	–	–	–	O ²⁾
Nobla upholstered		–	–	–	–	–	–	–	O
Nobla upholstered headboard, Nobla reha upholstered footboard		–	–	–	–	–	–	–	O
Decor									
Natural Beech ³⁾		✓	✓	✓	✓	–	–	–	✓
Royal Maple ³⁾		–	–	–	–	–	–	–	O
Havana Cherry ³⁾		–	–	–	–	–	–	–	O
Real Cherry		–	–	–	–	–	–	–	O ⁴⁾
Lindberg Oak		–	–	–	–	–	–	–	O
Bella Noce Choco ³⁾		–	–	–	–	–	–	–	–

Subject to change ✓ = Standard feature O = Available as an option – = Not available ¹⁾Only Econ with metal mattress base
²⁾Footboard also available as the lower Primus-Low version ³⁾Available for Relax senior bed ⁴⁾Only for Primus and Primus-Low foot/headboards

Model		Inovia	Westfalia IV	Westfalia-Klassik	Gigant	Allura II 100 120
Adjustments						
Height	Electric	✓	✓	✓	✓	✓
Backrest	Electric	✓	✓	✓	✓	✓
Thigh rest	Electric	✓	✓	✓	✓	✓
Locking feature	Selective	✓	✓	✓	✓	✓
Locking feature	Central	–	–	–	–	–
Wireless handset		–	–	–	–	–
24-volt drive system		✓	✓	✓	–	✓
Econ drive system		–	–	–	–	–
Height adjustment range						
from – to	cm	22-77	35-80	35-80	40-80	30-80
4-section mattress base						
Standard	cm	90 x 200 100 x 200	90 x 180 190 200	90 x 200	120 x 200	100 120 x 200
Special mattress base widths		–	100	–	–	–
Metal slats		✓	–	✓	–	–
Wire mesh mattress base		–	–	–	✓	✓
Sprung wooden slats		–	✓	–	–	–
Comfort mattress base		O	–	–	–	–
One-piece		✓	–	✓	✓	–
Multi-section		–	✓	–	–	✓
Adjustable mattress base tilt		✓	–	–	✓	✓
Sliding backrest		✓	–	✓	–	–
Ground clearance for lifting equipment		✓	✓	✓	✓	✓
Castor locking						
Central		–	–	✓	–	–
Per axle		✓	–	–	–	–
Single		–	✓	–	✓	✓
Castor diameter	mm	100	125	100	100	100
External dimensions and weight						
Length		220	212	212	222	222
Width		105 115	100	101	135	111 131
Net weight	kg	135 150	136	115	159	118 128
Safe working load						
	kg	225	185	200	350	250
Wooden surround						
Standard		–	✓ (90 cm only)	✓	–	–
Model-related headboard, closed surface		✓	–	–	✓	–
Model-related footboard, closed surface		O	–	–	✓	–
Model-related headboard with handle bar		O	–	–	–	✓
Model-related footboard with handle bar		✓	–	–	–	✓
Comodo with round posts		–	O	O	–	–
Primero with handle bar on head and footboard		–	O (90 cm only)	O	–	–
Prinzino with handle bar on head and footboard		–	–	–	–	–
Primus closed-surface headboard, footboard with handle bar		–	–	–	–	–
Nobla upholstered		–	–	–	–	–
Nobla upholstered headboard, Nobla reha upholstered footboard		–	–	–	–	–
Decor						
Natural Beech ²⁾		✓	✓	✓	✓	✓
Royal Maple ²⁾		–	–	–	–	–
Havana Cherry ²⁾		O	–	–	–	–
Real Cherry		–	–	–	–	–
Lindberg Oak		–	–	–	–	–
Bella Noce Choco ²⁾		–	–	–	–	–

Subject to change ✓ = Standard feature O = Available as an option – = Not available ¹⁾Only Econ with metal mattress base
²⁾Footboard also available as the lower Primus-Low version ³⁾Available for Relax senior bed ⁴⁾Only for Primus and Primus-Low foot/headboards

TECHNICAL DATA

Model	Dali	Dali wash	Dali low-econ	Dali low-entry	Lippe IV	Lippe IV washable	Lippe 120/140	Regia
Accessories								
Foam leather cover	O	O	O	O	–	–	–	O
Tray	O	O	O	O	–	–	–	O (DSG only)
Clamp-on extra-high safety side rail	O	O	O	–	–	–	–	–
Reading lamp	O	O	O	O	O	O	O	O
Bed extension 20 cm	O	–	O	–	–	–	–	O (27 cm)
Rastomat for lower leg rest	O ⁵⁾	✓	O	✓	✓	✓	✓	✓
Floor-mounted 24-volt power supply unit	O	O	–	O	O	O	O	–
Storage tray	O	O	O	O	–	–	–	O
Document holder	–	–	–	–	–	–	–	O
Transmitter holder for handset	O	O	O	O	O	O	–	O
Handset extension cable	–	–	–	–	O	O	–	O
Extension cushion	O	O	O	O	–	–	–	O (DSG only)
Linen holder	–	–	–	–	–	–	–	O
Wall deflection rollers	O	O	O	O	–	–	–	O
Wall spacers	–	–	–	–	–	–	–	O
Additional wooden surround (wooden panel)	O	–	O	O	–	–	–	–
Patient lifting pole including triangular handle	✓	✓	✓	✓	O	O	O	O
Mobilisation aid	O	–	O	O	O	O	O	–
Full-length safety side	✓	✓	✓	✓	–	–	–	O
Telescopic safety side	–	–	–	–	–	–	–	O
Upholstered side panels	–	–	–	–	–	–	–	O
Adaptable safety sides	O	–	O	–	O	–	–	–
Combinable safety side	–	–	–	–	–	–	–	–
Metal safety side	O	–	O	O	–	–	–	–
Lippe safety side	–	–	–	–	O	–	–	–
Transportation and storage aid	✓	✓	✓	✓	O	✓	–	–
Widening rail 5 cm	–	–	–	–	O	–	–	–
Tall mattress clamp handles (approx. 19 cm)	–	–	–	–	O	O	O	–
Softcovers	O	O	–	–	–	–	–	–
Carry bag for DSG bars	O	O	O	O	–	–	–	O


Model	Inovia	Westfalia IV	Westfalia-Klassik	Gigant	Allura II 100 120
Accessories					
Foam leather cover	O	O	O	O	O
Tray	O (DSG only)	O	O	–	O –
Clamp-on extra-high safety side rail	–	–	–	O (DSG only)	–
Reading lamp	O	O	O	O	O
Bed extension 20 cm	O (TSG: 27 cm)	O	O	O	O
Rastomat for lower leg rest	✓	✓	✓	✓	✓
Floor-mounted 24-volt power supply unit	–	O	O	–	–
Storage tray	O	O	O	O	O
Document holder	–	O	O	–	O
Transmitter holder for handset	–	O	O	–	–
Handset extension cable	–	O	O	–	–
Extension cushion	O (DSG only)	O	O	O (DSG only)	O
Linen holder	–	–	–	–	–
Wall deflection rollers	O	O	O	–	O
Wall spacers	–	–	–	–	–
Additional wooden surround (wooden panel)	–	–	–	–	O
Patient lifting pole including triangular handle	O	✓	✓	✓	✓
Mobilisation aid	–	O	–	–	O
Full-length safety side	✓	✓	✓	✓	✓
Telescopic safety side	O	–	–	–	–
Upholstered side panels	–	–	–	–	–
Adaptable safety sides	–	O	–	–	–
Combinable safety side	–	–	–	O	–
Metal safety side	–	–	–	–	–
Lippe safety side	–	–	–	–	–
Transportation and storage aid	–	–	–	–	✓
Widening rail 5 cm	–	–	–	–	–
Raised mattress clamp handles (approx. 19 cm)	–	–	–	–	–
Softcovers	–	–	–	–	–
Carry bag for DSG bars	–	–	–	–	–

Dali range



Special features

- Wireless handset for added comfort
- Apps provide digital support
- Full-length safety sides with Easy Click system
- Machine washable version
- 24-volt system dispenses with leakage current testing

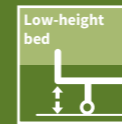

Lippe range/ Relax



Special features

- Quickly assembled by just one person
- Small supporting frame for effortless installation in almost all bed frames
- Height adjustment from 29 - 74 cm
- Machine washable version




Regia/Inovia



Special features

- Low-height bed for fall prevention
- Comfort and elegance for private customers
- Split safety sides can be selected
- Choice of homelike decors, fabrics and faux leather

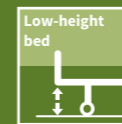

Westfalia range



Special features

- Care beds with stable scissor-type frame
- Attractive head and footboard designs
- Selectively lockable handset
- Cost-saving 24-volt drive system



Heavy-duty beds



Special features

- High safe working load, up to 350 kg
- Mattress base widths up to 120 cm
- Better quality of life for heavier residents



Vertica homecare



Special features

- Innovative mobilisation bed
- Genuine sitting position
- Ergonomic position for standing up
- Ideal for care at home




Bedside cabinets and side tables



Special features

- Storage space for personal belongings and care utensils
- Ideal for drinks and meals taken in bed
- Convenient bookrests
- Good manoeuvrability for flexible use

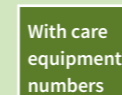

Transport chair/ mattresses/ furniture



Special features

- Customised furnishings for care rooms
- Easily accessible storage space
- Save and comfortable transfer in the Ravello transport chair

Dali range



Dali range





Quick and digital


Dali

ADVANTAGES

The Dali care bed combines comfort, safety and easy operation at an excellent price/performance ratio. The Easy Click system allows simple and effortless assembly and dismantling. A wireless handset and helpful apps open the door to the digital future of care at home.

- Safe working load: 185 kg
- Maximum weight of resident: 155 kg
- Height adjustment range: 40 to 80 cm
- 24-volt drive system
- No leakage current measurements during expected service life
- Reverse-Trendelenburg position (optionally without)
- Safety side installation with Easy Click system requires no tools
- Wireless handset
- Dali-Lock app for managing the adjustment functions

External dimensions/height adjustment range:

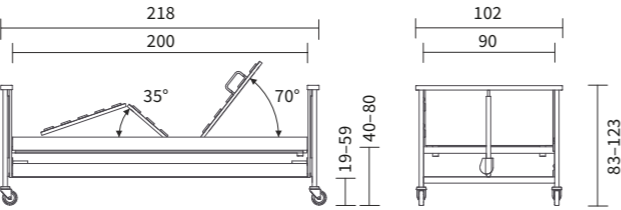


Special features:



Wireless handset for greater freedom of movement

Safety side protective height of 37 cm

185 kg safe working load, 4-section mattress base

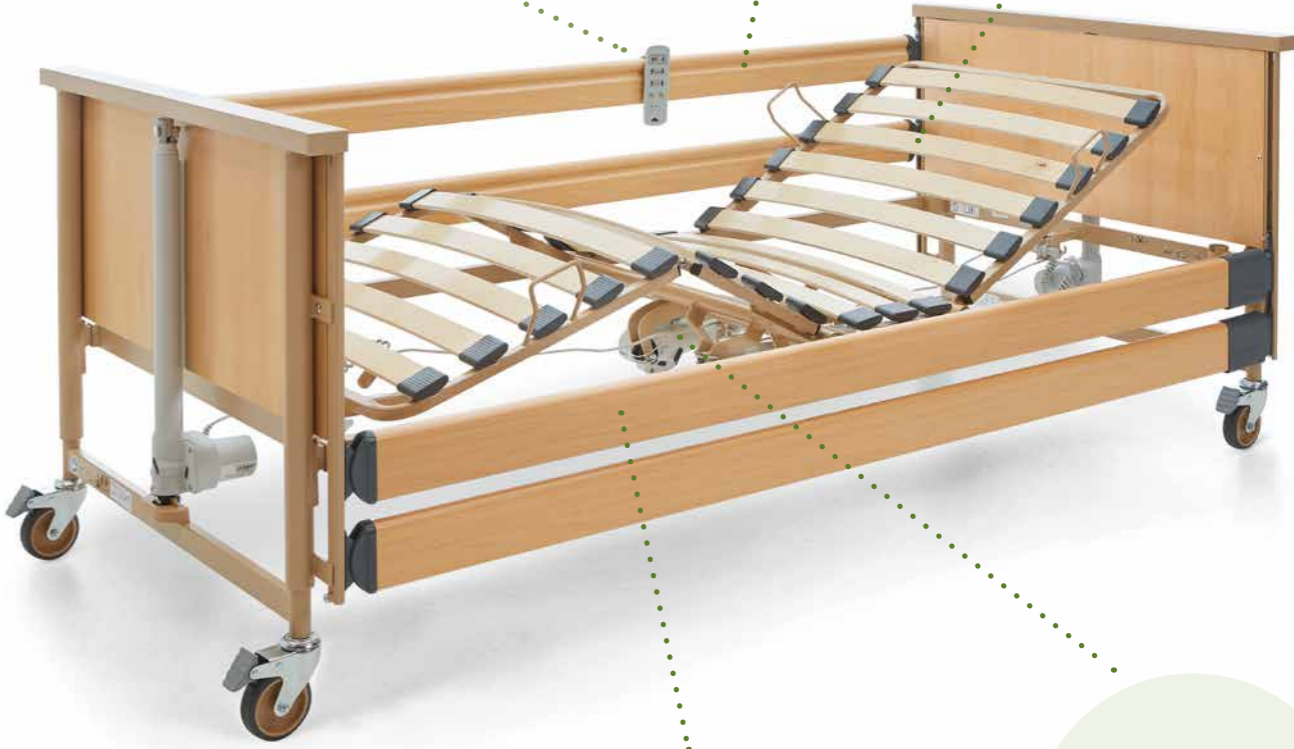

Castors with easy braking

Assembly of safety sides with Easy Click system requires no tools

Cost-saving 24-volt drive system

Further information:



Product page

Assembly video

Info video

The advantages at a glance

DESIGN DETAILS



Reverse-Trendelenburg position.



Safety sides with Easy Click system.



Compact delivery and storage on storage aid/transportation aid.



Dali-Lock app for managing all adjustment functions.

ACCESSORIES



Bed extension including safety side bars (also for retrofitting).



The clamp-on safety side bar raises the side protection height by 18 cm.



Soft foam leather covers protect the resident when in contact with the safety sides and when leaving the bed.



With the help of the optional ratchet mechanism, the lower leg rest can be manually raised.



The adaptable safety side for Dali beds offers good protection and leaves enough space to get in and out of bed.



The metal safety side is particularly robust and easy to clean.



The tray is simply placed on top of the raised safety side bars.



The storage tray for glasses, mobiles or the wireless handset can be attached at various positions.



Wall deflection rollers protect the bed and the walls from damage.



The extension cushion can be fitted onto the full-length safety side to make it easier for residents to sit on the edge of the bed.



The reading lamp can be fitted to the patient lifting pole or into the connection sleeve.



Elegant Softcovers for head and footboards and safety sides.



Homelike wooden surrounds for the head and footboards.



With the floor-mounted power supply unit, the 24-volt system can also be used with various international mains cables.



The Dali low-entry-econ model has a safe Econ 230-volt drive system.



The Econ drive system includes a corded handset with selectively lockable functions.

EQUIPMENT AND OPTIONS

Equipment/Option	Dali	Material no.
Design features		
Dali with sprung wooden slat mattress base and wooden safety side		276349
Dali with metal slatted mattress base and wooden safety side		276354
Dali with sprung wooden slat mattress base, wooden safety side and wood surround head/footboards		276358
Dali with metal slatted mattress base, wooden safety side and wood surround head/footboards		276359
Dali econ with sprung wooden slat mattress base and wooden safety side		283947
Dali econ with metal slatted mattress base and wooden safety side		283948
Dali econ with metal slatted mattress base and wooden safety side 90 x 180 cm		281206
Technical data		
Net weight	97 kg	
Overall dimensions	102 x 218 cm	
Height adjustment range	40–80 cm	
Safe working load	185 kg	
Maximum resident weight	145 – 155 kg	
24-volt system	✓	
Prepared for use of Dali-Lock app	✓ (24-volt only)	
Econ 230-volt drive system	0	
Standard features		
Frame colour	Topaz	
Decor	Natural Beech	
Mattress base size	90 x 200 cm	
Sprung wooden slat mattress base (mattress base model)	✓	
Full-length safety side (DSG)	✓	
Type of handset (selectively lockable/wireless)	24-volt ✓/✓	Econ ✓/–
Mechanical CPR release	✓	
Castor diameter	100 mm	
Castor locking	Single	
Patient lifting pole with triangular grab handle	✓	
Special equipment and accessories		
Full-length metal safety side	1 bar	271368
Adaptable safety sides	Right	250138
	Left	250139
Clamp-on safety side bar	0	191367
Wood surrounds as chassis cladding	2 pieces per set	282157
Mobilisation aid	0	175402
Bed extension for sprung wooden slat mattress base + safety side bars	0	1 x 266002
		4 x 271363
Bed extension for metal slatted mattress base + safety side bars	0	1 x 266003
		4 x 271363
Mattress piece for bed extension	0	166715
Wall deflection roller	0	206752
Foam leather cover for full-length safety sides (DSG)	1 side	225921
Reading lamp for easy insertion	0	264877
Reading lamp for fixing to the patient lifting pole	0	1 x 191879
		1 x 174231
Under bed light		
Tray for full-length safety sides (DSG)	0	174296
Extension cushion	0	258395
Transmitter holder for handset	0	199137
Carry bag for safety side bars	0	273831
Softcover	Vanilla	259579
	Bordeaux	259579
	Orange	259579
	Petrol	259579
	Dark blue	259579
	Lavender	259579
Ratchet mechanism, electrical adjustment of lower leg rest for sprung wooden slat mattress base	0	197091
Storage tray	White	266970
	Yellow	266971
	Red	266619
Floor-mounted 24-volt power supply unit	0	254276

✓ = Standard feature 0 = Available as an option – = Not available



The bed for optimal hygiene


Dali wash

ADVANTAGES

The Dali wash care bed was designed to allow cleaning in automated washing systems. With its high protection category of IPX6, the bed can be reprocessed without uninstalling the drive system. The bed is protected from corrosion by its cathodic dip coating (CDC) and holes that allow residual fluid to drain away.

- Safe working load: 185 kg
- Maximum weight of resident: 155 kg
- Height adjustment range: 40 to 80 cm
- 24-volt drive system
- No leakage current measurements during expected service life
- Reverse-Trendelenburg position (optionally without)
- Machine washable for optimal hygiene
- Safety side installation with Easy Click system requires no tools
- Wireless handset
- Dali-Lock app for managing the adjustment functions

External dimensions/height adjustment range:

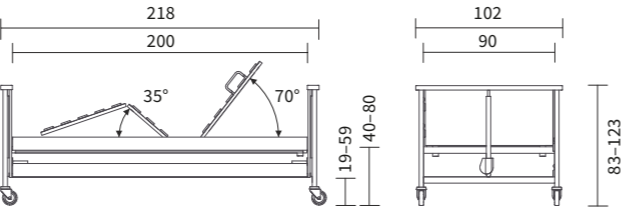

Special features:



Further information:



Product page




Assembly video



Info video



Machine washable

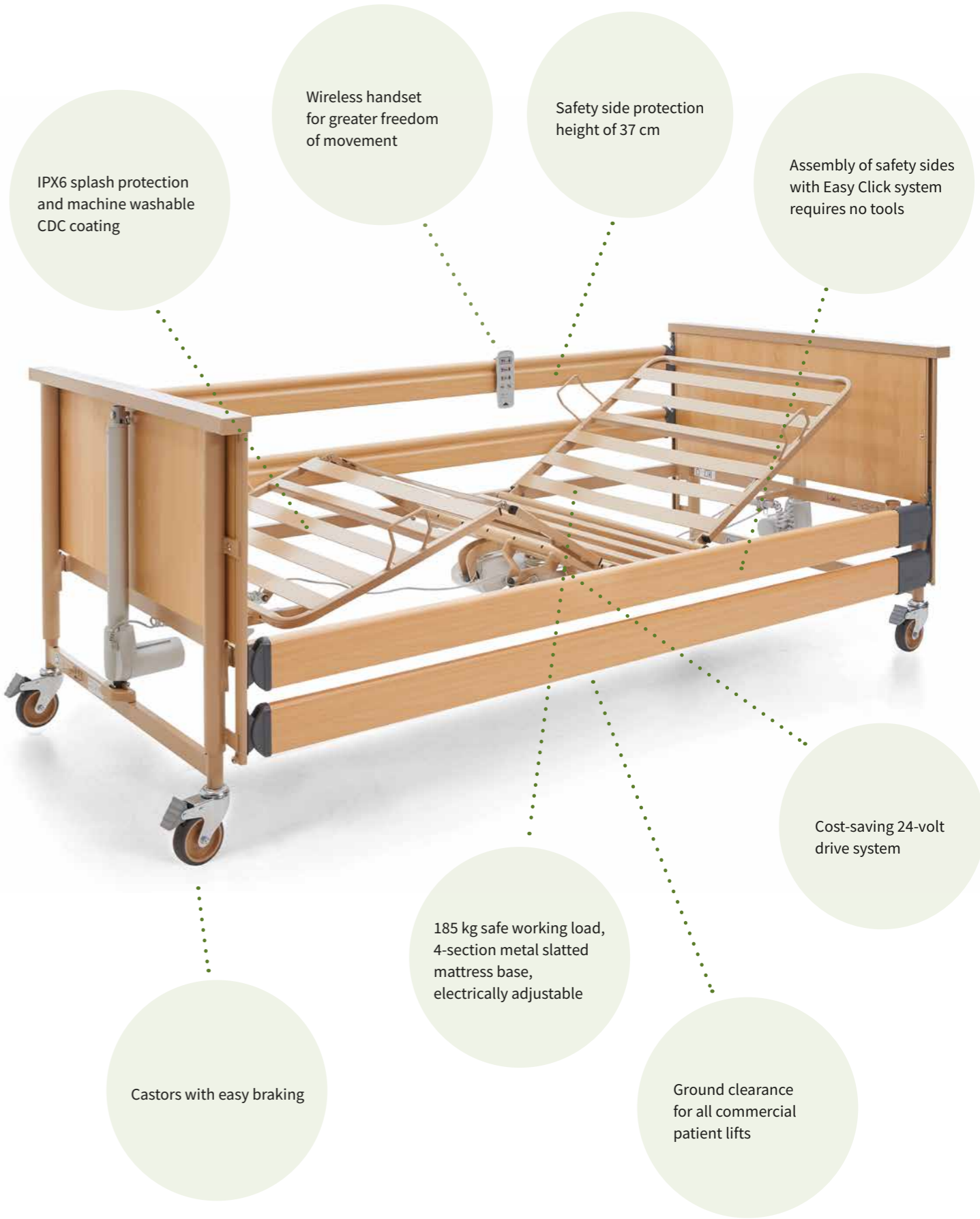

The advantages at a glance

DESIGN DETAILS



The CDC coating protects the bed from corrosion, inside and out.



Drilled holes in the metal profiles allow residual fluid to drain away.



Plastic plugs with drilled holes allow residual moisture to escape.



The sealed surfaces are ideal for safe and thorough cleaning.



The 4-section mattress base has hygienic metal slats.

ACCESSORIES



The clamp-on safety side bar raises the side protection height by 18 cm.



Soft foam leather covers protect the resident when in contact with the safety sides and when leaving the bed.



The bars of the safety sides can be transported and stored even more easily with the carry bag.



The tray is simply placed on top of the raised safety side bars.



The storage tray for glasses, mobiles or the wireless handset can be attached at various positions.



The extension cushion can be fitted onto the full-length safety side to make it easier for residents to sit on the edge of the bed.



The reading lamp can be fitted to the patient lifting pole or into the connection sleeve.



Elegant Softcovers for head and footboards and safety sides.



With the floor-mounted power supply unit, the 24-volt system can also be used with various international mains cables.

EQUIPMENT AND OPTIONS

Equipment/Option	Dali wash	Material no.
Design features		
Dali wash with metal slatted mattress base and wooden safety side		278060
Technical data		
Net weight	97 kg	
Overall dimensions	102 x 218 cm	
Height adjustment range	40 – 80 cm	
Safe working load	185 kg	
Maximum resident weight	155 kg	
24-volt system	✓	
Prepared for use of Dali-Lock app	✓	
Standard features		
Frame colour	Topaz	
Decor	Natural Beech	
Mattress base size	90 x 200 cm	
Metal slatted mattress base (mattress base model)	✓	
Safety side	Full-length safety side (DSG)	
Type of handset (selectively lockable/wireless)	✓	
Mechanical CPR release	✓	
Castor diameter	100 mm	
Castor locking	Single	
Ratchet mechanism, electrical adjustment of lower leg rest for sprung wooden slat mattress base	✓	
Patient lifting pole with triangular grab handle	✓	
Special equipment and accessories		
Clamp-on safety side bar	O	191367
Foam leather cover for full-length safety sides (DSG)	1 side	225921
Reading lamp for easy insertion	O	264877
Reading lamp for fixing to the patient lifting pole	O	1 x 191879 1 x 174231
Tray for full-length safety sides (DSG)	O	174296
Extension cushion	O	258395
Transmitter holder for handset	O	199137
Carry bag for safety side bars	O	273831
Softcover	Vanilla	259579
	Bordeaux	259579
	Orange	259579
	Petrol	259579
	Dark blue	259579
	Lavender	259579
Storage tray	White	266970
	Yellow	266971
	Red	266619
Floor-mounted 24-volt power supply unit	O	254276

✓ = Standard feature O = Available as an option



The budget-oriented low-height bed


Dali low-econ

ADVANTAGES


The Dali low-econ combines many strengths at an excellent price/performance ratio. With a height adjustment range from 32 to 72 cm, its lowest position protects residents from falls. It also allows taller people to maintain a back-friendly, upright posture when working at the bed.

- Safe working load: 185 kg
- Maximum weight of resident: 155 kg
- Height adjustment range: 32 to 72 cm
- Econ drive system
- Reverse-Trendelenburg position
- Safety side installation with Easy Click system requires no tools
- Selective lockable, cabled handset

External dimensions/height adjustment range:

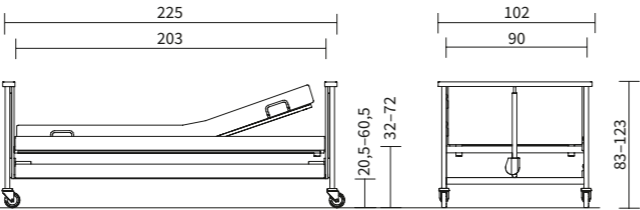

Special features:

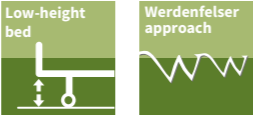

Further information:



DESIGN DETAILS



Mattress base with sprung wooden slats.



Compact delivery and storage on storage aid/transportation aid.



Safe and reliable Econ 230-volt drive system.



Selectively lockable handset.

ACCESSORIES



Bed extension including safety side bars (also for retrofitting).



The clamp-on safety side bar raises the side protection height by 18 cm.



Soft foam leather covers protect the resident when in contact with the safety sides and when leaving the bed.



The bars of the safety sides can be transported and stored even more easily with the carry bag.



With the help of the optional ratchet mechanism, the lower leg rest can be manually raised.



The adaptable safety side for Dali beds offers good protection and leaves enough space to get in and out of bed.



The metal safety side is particularly robust and easy to clean.



The tray is simply placed on top of the raised safety side bars.



The storage tray for glasses, mobiles or the wireless handset can be attached at various positions.



Wall deflection rollers protect the bed and the walls from damage.



The extension cushion can be fitted onto the full-length safety side to make it easier for residents to sit on the edge of the bed.



The reading lamp can be fitted to the patient lifting pole or into the connection sleeve.




Homelike wooden surrounds for the head and footboards.

EQUIPMENT AND OPTIONS

Equipment/Option	Dali low-econ	Material no.
Design features		
Dali low-econ with sprung wooden slat mattress base and wooden safety side		283952
Dali low-econ with sprung wooden slat mattress base and wooden safety side with wood surround head/footboards		283951
Technical data		
Net weight	97 kg	
Overall dimensions	102 x 225 cm	
Height adjustment range	32–72 cm	
Safe working load	185 kg	
Maximum resident weight	145 – 155 kg	
Econ 230-volt drive system	✓	
Standard features		
Frame colour	Topaz	
Decor	Natural Beech	
Mattress base size	90 x 200 cm	
Sprung wooden slat mattress base (mattress base model)	✓	
Full-length safety side (DSG)	✓	
Type of handset (selectively lockable)	✓	
Mechanical CPR release	✓	
Castor diameter	100 mm	
Castor locking	Single	
Patient lifting pole with triangular grab handle	✓	
Special equipment and accessories		
Full-length metal safety side	1 bar	271368
Adaptable safety sides	Right	250138
	Left	250139
Clamp-on safety side bar	O	191367
Wood surrounds as chassis cladding	2 pieces per set	282159
Mobilisation aid	O	175402
Bed extension for sprung wooden slat mattress base + safety side bars	O	1 x 266002
		4 x 271363
Mattress piece for bed extension	O	166715
Wall deflection roller	O	206752
Foam leather cover for full-length safety sides (DSG)	1 side	225921
Reading lamp for easy insertion	O	264877
Reading lamp for fixing to the patient lifting pole	O	1 x 191879
		1 x 174231
Tray for full-length safety sides (DSG)	O	174296
Extension cushion	O	258395
Transmitter holder for handset	O	199137
Carry bag for safety side bars	O	273831
Ratchet mechanism, electrical adjustment of lower leg rest	O	197091
Storage tray	White	266970
	Yellow	266971
	Red	266619

✓ = Standard feature O = Available as an option

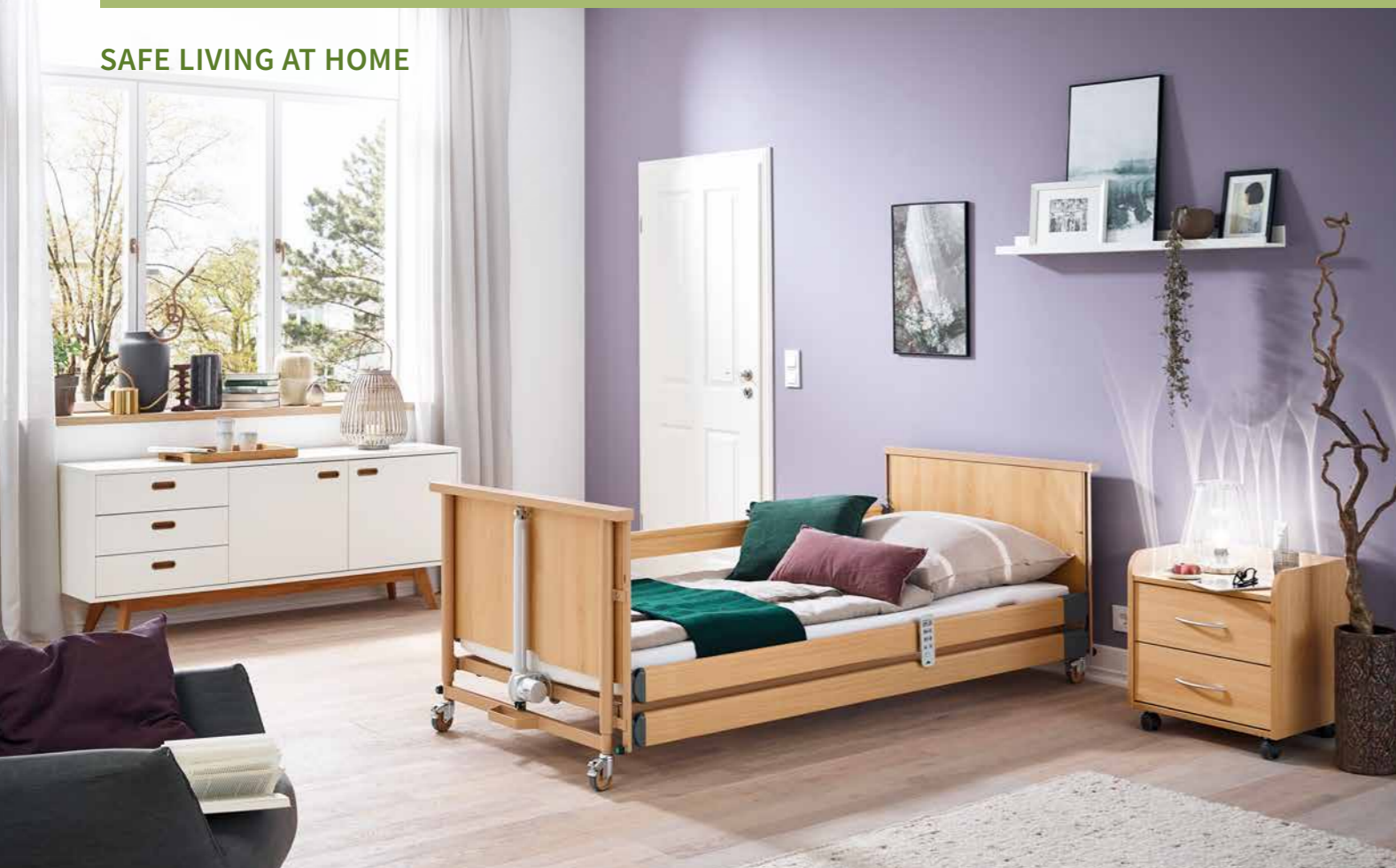

The low-height bed


Dali low-entry

ADVANTAGES

The low-height bed Dali low-entry prevents falls without resorting to deprivation of liberty safeguards. With its height adjustment range of 23 to 63 cm, it is close to the floor when at its lowest position. The risk of injury is then much lower in the event of a fall. This makes it easier to decide not to use safety sides.

- Safe working load: 175 kg
- Maximum weight of resident: 145 kg
- Height adjustment range: 23 to 63 cm
- 24-volt drive system
- No leakage current measurements during expected service life
- Reverse-Trendelenburg position (optionally without)
- Safety side installation with Easy Click system requires no tools
- Wireless handset
- Dali-Lock app for managing the adjustment functions

External dimensions/height adjustment range:

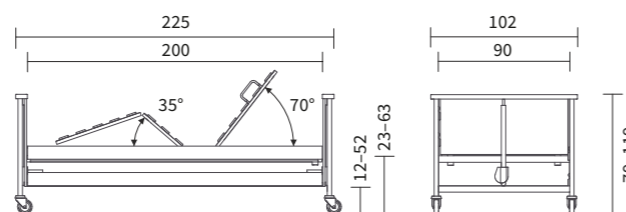

Special features:



Low-level entry thanks to height adjustment from 23 to 63 cm

Wireless handset for greater freedom of movement

175 kg safe working load
4-section metal slatted mattress base, electrically adjustable

Assembly of safety sides with Easy Click system requires no tools



Ratchet mechanism for adjusting the lower leg rest

Reverse-Trendelenburg position can be set

Cost-saving 24-volt drive system

Further information:



DESIGN DETAILS



Reverse-Trendelenburg position.



Electronic (SwitchMode) power supply unit integrated immediately downstream of the mains plug.



Power socket for connecting the electronic power supply unit (SwitchMode).



Compact delivery and storage on storage aid/transportation aid.

ACCESSORIES



Soft foam leather covers protect the resident when in contact with the safety sides and when leaving the bed.



The bars of the safety sides can be transported and stored even more easily with the carry bag.



The tray is simply placed on top of the raised safety side bars.



The storage tray for glasses, mobiles or the wireless handset can be attached at various positions.



Wall deflection rollers protect the bed and the walls from damage.



The extension cushion can be fitted onto the full-length safety side to make it easier for residents to sit on the edge of the bed.



Reading lamps.



Homelike wooden surrounds for the head and footboards.



The metal safety side is particularly robust and easy to clean.



With the floor-mounted power supply unit, the 24-volt system can also be used with various international mains cables.



The Dali low-entry-econ model has a safe Econ 230-volt drive system.



The Econ drive system includes a corded handset with selectively lockable functions.

EQUIPMENT AND OPTIONS

Equipment/Option	Dali low-entry	Material no.
Design features		
Dali low-entry with metal slatted mattress base and wooden safety side		276356
Dali low-entry-econ with metal slatted mattress base and wooden safety side		283954
Technical data		
Net weight	97 kg	
Overall dimensions	102 x 225 cm	
Height adjustment range	23 – 63 cm	
Safe working load	175 kg	
Maximum resident weight	145 kg	
24-volt system	✓	
Prepared for use of Dali-Lock app	✓ (24-volt only)	
Econ 230-volt drive system	0	
Standard features		
Frame colour	Topaz	
Decor	Natural Beech	
Mattress base size	90 x 200 cm	
Metal slatted mattress base (mattress base model)	✓	
Safety side	Full-length safety side (DSG)	
Type of handset (selectively lockable/wireless)	24-volt ✓/✓	Econ ✓/–
Mechanical CPR release		✓
Castor diameter	75 mm	
Castor locking	Single	
Ratchet mechanism, electrical adjustment of lower leg rest for sprung wooden slat mattress base		✓
Patient lifting pole with triangular grab handle		✓
Special equipment and accessories		
Full-length metal safety side	1 bar	271368
Wood surrounds as chassis cladding	2 pieces per set	282159
Wall deflection roller	1 piece	206752
Foam leather cover for full-length safety sides (DSG)	1 side	225921
Reading lamp for easy insertion	0	264877
Reading lamp for fixing to the patient lifting pole	0	1 x 191879 1 x 174231
Tray for full-length safety sides (DSG)	0	174296
Extension cushion	0	258395
Transmitter holder for handset	0	199137
Carry bag for safety side bars	0	273831
Storage tray	White	266970
	Yellow	266971
	Red	266619
	0	254276
Floor-mounted 24-volt power supply unit	0	254276


✓ = Standard feature 0 = Available as an option – = Not available

Lippe range/Relax



Lippe range/Relax



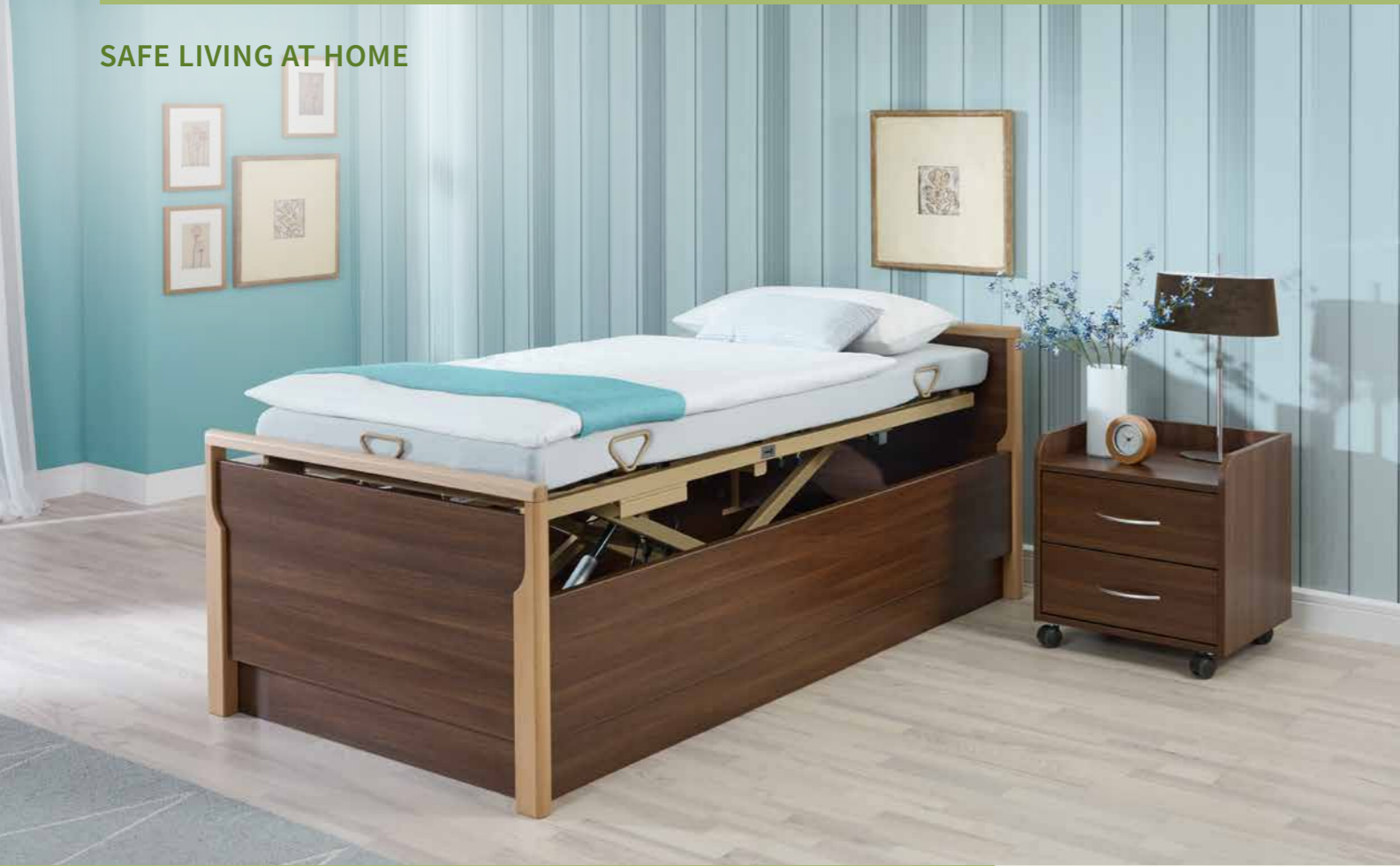

The height-adjustable bed insert


Lippe IV

ADVANTAGES



The Lippe IV bed insert can be installed in an existing bed frame and offers all the advantages of an electrically adjustable care bed. The large height adjustment range makes getting out of bed easy and allows carers to work at an ergonomically ideal height.

- Safe working load: 200 kg
- Maximum weight of resident: 170 kg
- Height adjustment range: 29 to 74 cm
- 24-volt drive system
- No leakage current measurements during expected service life
- Mattress base dimensions of 80, 90 and 100 cm in width and 190 and 200 cm in length
- Space saving supporting frame measuring only 69 x 113 cm
- Secure and stable thanks to feet with floor-level adjustment screws

External dimensions/height adjustment range:

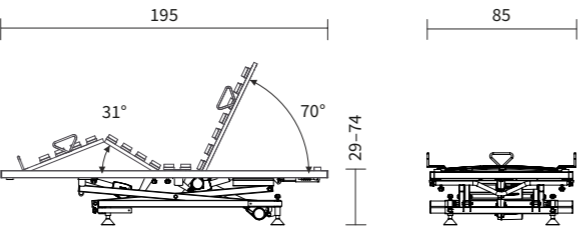

Special features:

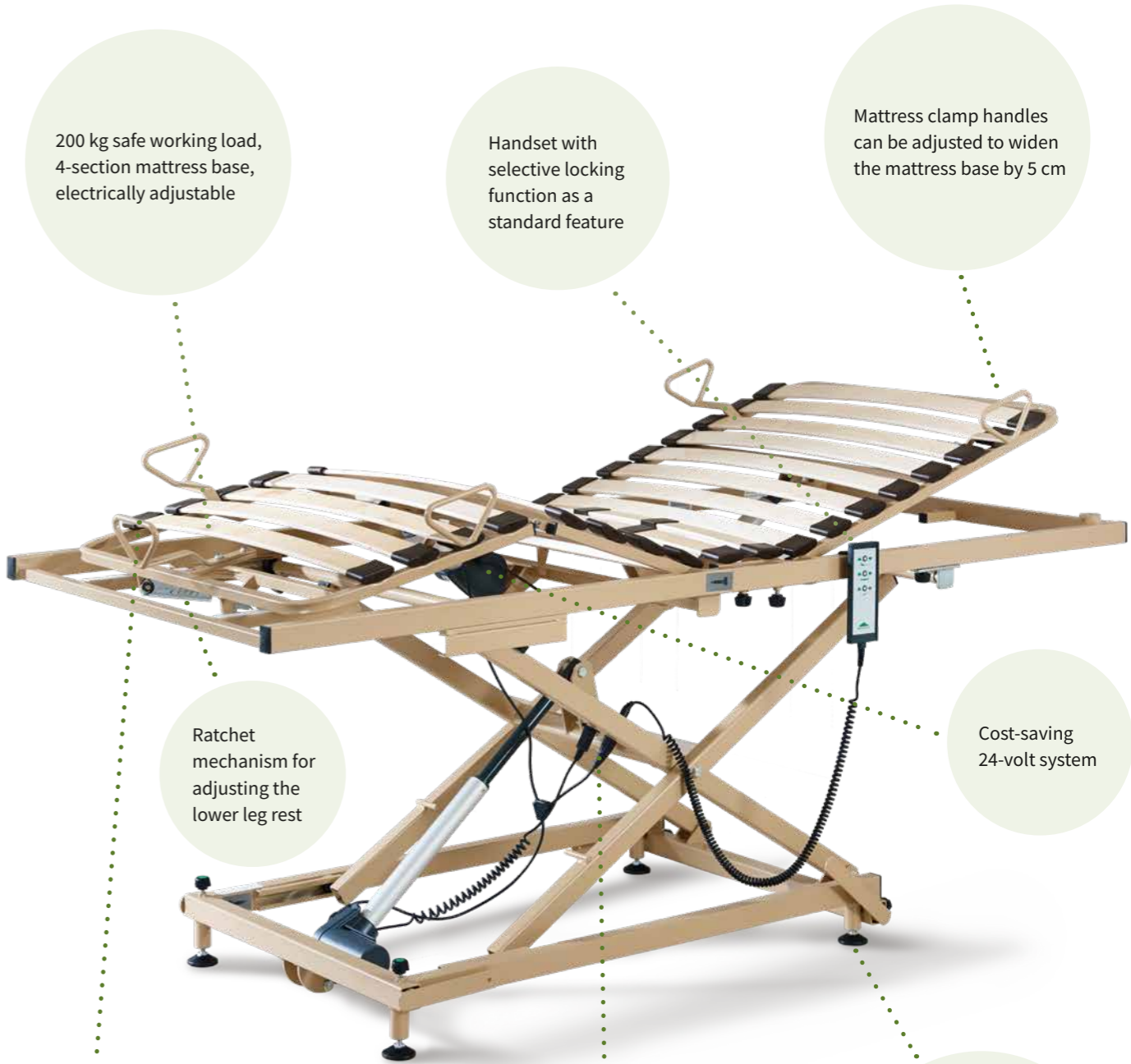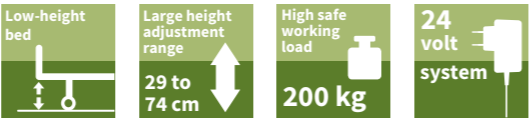

Further information:



Product page



Assembly video

The advantages at a glance

DESIGN DETAILS



Removable mattress clamp handles for lengthening and widening by 5 cm in each case.



Electronic (SwitchMode) power supply unit integrated immediately downstream of the mains plug.



Selectively lockable handset.



Feet with floor-level adjustment screws and fixing lugs.



Power socket for connecting the electronic power supply unit (SwitchMode).

ACCESSORIES



Storage and transport aid with steerable castors.



Reading lamps.



Patient lifting pole.



Lippe safety side (also for retrofitting).



Mobilisation aid.



Adjustable feet for raising the lowest position (+ 7 cm).



Widening rail, 5 cm.



Adaptable safety sides.



Taller mattress clamp handles (approx. 19 cm).



Flexible transmitter holder for the handset.



Extension cable for the handset (1.40 m).



With the floor-mounted power supply unit, the 24-volt system can also be used with various international mains cables.

EQUIPMENT AND OPTIONS

Equipment/Option	Lippe IV	Material no.
Design features		
Lippe IV 90 x 190 cm		278476
Lippe IV 90 x 200 cm		278478
Lippe IV 100 x 190 cm		275091
Lippe IV 100 x 200 cm		275092
Lippe IV 80 x 190 cm		275087
Lippe IV 80 x 200 cm		275088
Technical data		
Net weight	80 kg	
Overall dimensions	85 x 185 cm	
	85 x 195 cm	
	95 x 185 cm	
	95 x 195 cm	
	75 x 185 cm	
	75 x 195 cm	
Height adjustment range	29 -74 cm	
Safe working load	200 kg	
Maximum resident weight	170 kg	
24-volt system	✓	
Standard features		
Frame colour	Topaz	
Mattress base size	90 x 200 cm	
Sprung wooden slat mattress base (mattress base model)	✓	
Type of handset (selectively lockable/wireless)	Selective	
Mechanical CPR release	✓	
Ratchet mechanism, electrical adjustment of lower leg rest	✓	
Special equipment and accessories		
Adaptable safety sides	Right	250138
	Left	250139
Lippe safety side		180457
Adjustable feet	1 piece	180431
Widening rail for mattress base	O	183505
Patient lifting pole with triangular grab handle	O	146602
Mobilisation aid	O	175402
Reading lamp for easy insertion	O	264877
Reading lamp for fixing to the patient lifting pole	O	1 x 191879
		1 x 174231
Transmitter holder for handset	O	199137
Handset extension cable	O	178986
Transportation and storage aid	O	Right: 183724 +
		Left: 183725
Floor-mounted 24-volt power supply unit	O	254276

✓ = Standard feature O = Available as an option



The bed insert for optimal hygiene


Lippe IV washable

ADVANTAGES

The Lippe IV washable was designed to allow cleaning in automated washing systems. With its high protection category of IPX6, the bed insert can be reprocessed without uninstalling the drive system. It is protected from corrosion by its cathodic dip coating (CDC) and holes that allow residual fluid to drain away.

- Safe working load: 200 kg
- Maximum weight of resident: 170 kg
- Height adjustment range: 29 to 74 cm
- 24-volt drive system
- No leakage current measurements during expected service life
- Machine washable for optimal hygiene
- Mattress base size 90 x 190 cm and 90 x 200 cm
- Selectively lockable handset
- Secure and stable thanks to feet with floor-level adjustment screws

External dimensions/height adjustment range:

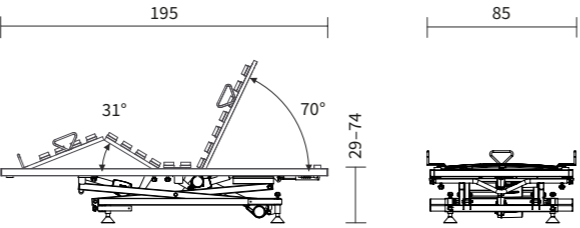

Special features:


Machine washable

Low-height bed

Large height adjustment range
29 to 74 cm

High safe working load
200 kg

24 volt system

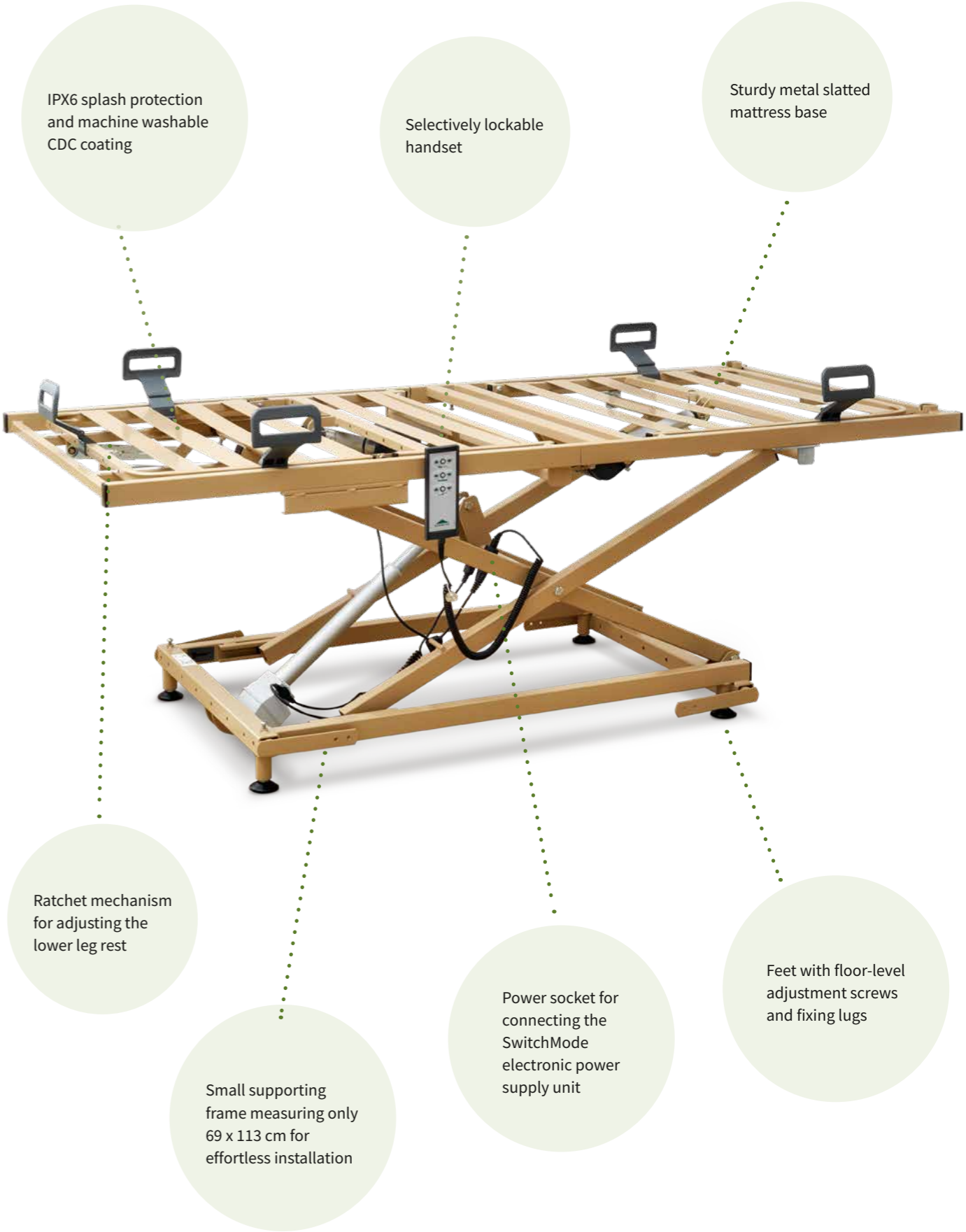

Further information:



The advantages at a glance

DESIGN DETAILS



Removable mattress clamp handles for lengthening and widening by 5 cm in each case.



Electronic (SwitchMode) power supply unit integrated immediately downstream of the mains plug.



Selectively lockable handset.



Feet with floor-level adjustment screws and fixing lugs.



Power socket for connecting the electronic power supply unit (SwitchMode).



Storage and transport aid with steerable castors.

ACCESSORIES



Reading lamps.



Patient lifting pole.



Mobilisation aid.



Adjustable feet for raising the lowest position (+ 7 cm).



Taller mattress clamp handles (approx. 19 cm).



Flexible transmitter holder for the handset.



Extension cable for the handset (1.40 m).




With the floor-mounted power supply unit, the 24-volt system can also be used with various international mains cables.

EQUIPMENT AND OPTIONS

Equipment/Option	Lippe IV washable	Material no.
Design features		
Lippe IV washable 90 x 190 cm		275097
Lippe IV washable 90 x 200 cm		275098
Technical data		
Net weight	80 kg	
Overall dimensions	85 x 185 cm	
	85 x 195 cm	
Height adjustment range	29 -74 cm	
Safe working load	200 kg	
Maximum resident weight	170 kg	
24-volt system	✓	
Standard features		
Frame colour	Topaz	
Mattress base size	90 x 190 cm	
	90 x 200 cm	
Metal slatted mattress base (mattress base model)	✓	
Type of handset (selectively lockable/wireless)	Selective	
Mechanical CPR release	✓	
Ratchet mechanism, electrical adjustment of lower leg rest	✓	
Special equipment and accessories		
Adjustable feet	1 piece	180431
Patient lifting pole with triangular grab handle	O	146602
Mobilisation aid	O	175402
Reading lamp for easy insertion	O	264877
Reading lamp for fixing to the patient lifting pole	O	1 x 191879
		1 x 174231
Transmitter holder for handset	O	199137
Handset extension cable	O	178986
Floor-mounted 24-volt power supply unit	O	254276

✓ = Standard feature O = Available as an option

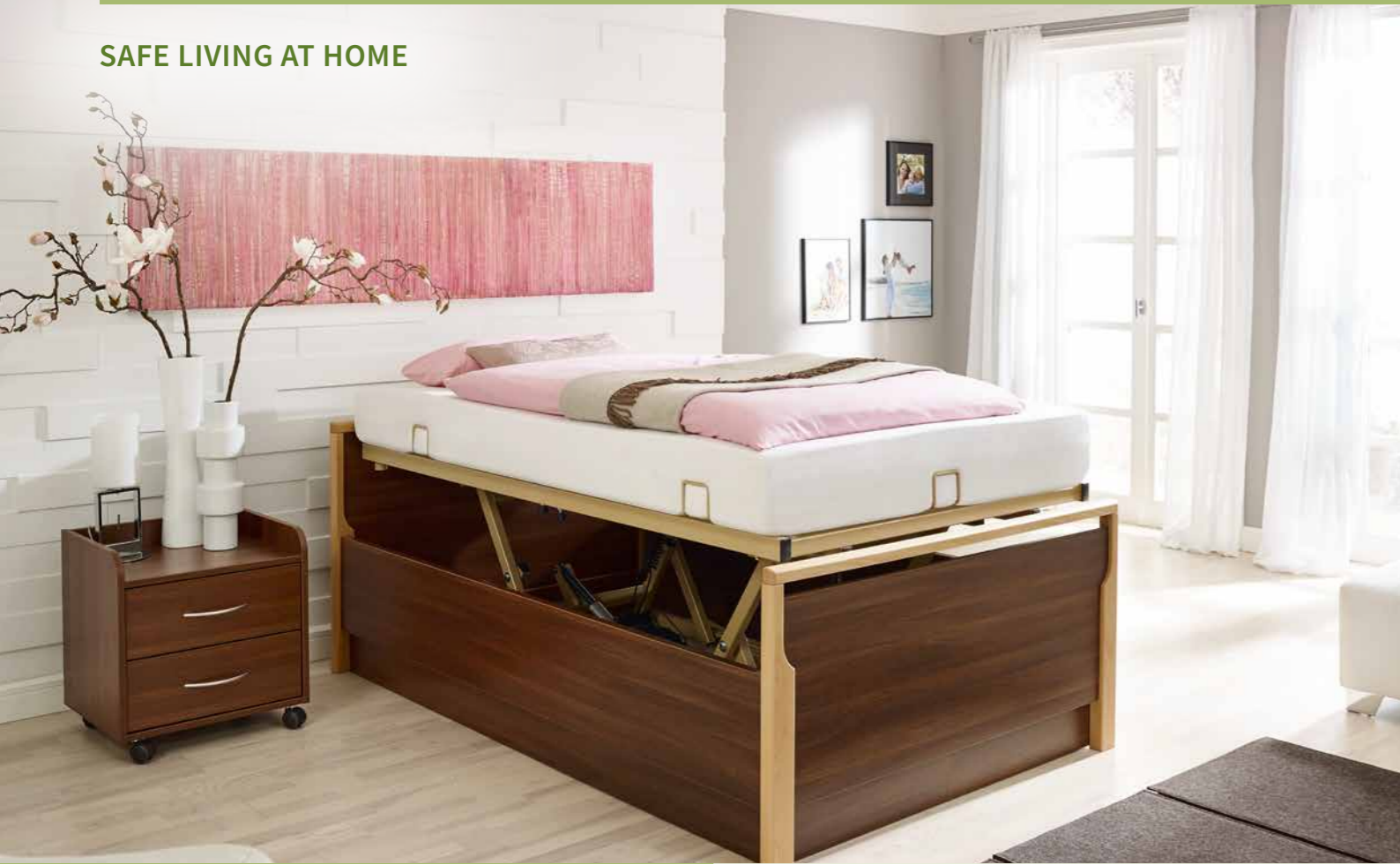

Plenty of space and comfort
in any circumstances


Lippe 120/140

ADVANTAGES


The Lippe 120 and Lippe 140 bed inserts with mattress base widths of 120 and 140 cm offer plenty of space and comfort. With their high safe working load of 250 kg, they are also a practical alternative to a heavy-duty bed. Depending on the type of mattress and accessories used, these models can accommodate residents who weigh more than 200 kg.

- Safe working load: 250 kg
- Maximum weight of resident: 210 kg
- Height adjustment range: 29 to 74 cm
- 24-volt drive system
- No leakage current measurements during expected service life
- Selectively lockable handset
- Secure and stable thanks to feet with floor-level adjustment screws

External dimensions/height adjustment range:

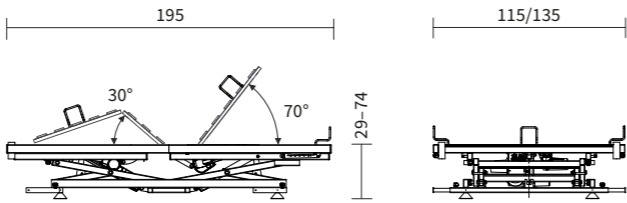

Special features:

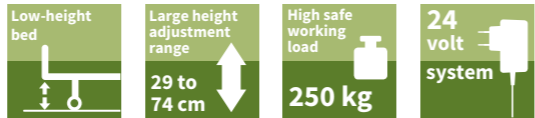

Further information:



Product page



Assembly video

The advantages at a glance

DESIGN DETAILS



Feet with floor-level adjustment screws and fixing lugs.



Height adjustment with ball bearings.



Electronic (SwitchMode) power supply unit integrated immediately downstream of the mains plug.



Power socket for connecting the electronic power supply unit (SwitchMode).



The low entry height of 29 cm makes getting in and out of bed easier.

ACCESSORIES



Reading lamps.



Mobilisation aid.



Adjustable feet for raising the lowest position (+ 7 cm).



Bed extension (also for retrofitting).



Patient lifting pole.



Taller mattress clamp handles (approx. 19 cm).




With the floor-mounted power supply unit, the 24-volt system can also be used with various international mains cables.

EQUIPMENT AND OPTIONS

Equipment/Option	Lippe 120/140	Material no.
Design features		
Lippe 120 x 200 cm		275093
Lippe 120 x 220 cm		275095
Lippe 140 x 200 cm		275094
Lippe 140 x 220 cm		275096
Technical data		
Net weight	129/138 kg	
Overall dimensions	115 x 195 cm 135 x 195 cm	
Height adjustment range	29 -74 cm	
Safe working load	250 kg	
Maximum resident weight	210 kg	
24-volt system	✓	
Standard features		
Frame colour	Topaz	
Mattress base size	120 x 200 cm 140 x 200 cm	
Metal slatted mattress base (mattress base model)	✓	
Type of handset (selectively lockable/wireless)	Selective	
Mechanical CPR release	✓	
Special equipment and accessories		
Mobilisation aid	O	244029
Patient lifting pole with triangular grab handle	O	146602
20-cm bed length extension for 120-cm wide metal slatted mattress base	O	240285
20-cm bed length extension for 140-cm wide metal slatted mattress base	O	240286
Mattress piece for bed extension, 120 cm wide	O	241942
Mattress piece for bed extension, 140 cm wide	O	244859
Reading lamp for easy insertion	O	264877
Reading lamp for fixing to the patient lifting pole	O	1 x 191879 1 x 174231
Handset extension cable	O	238047
Floor-mounted 24-volt power supply unit	O	254276

✓ = Standard feature O = Available as an option

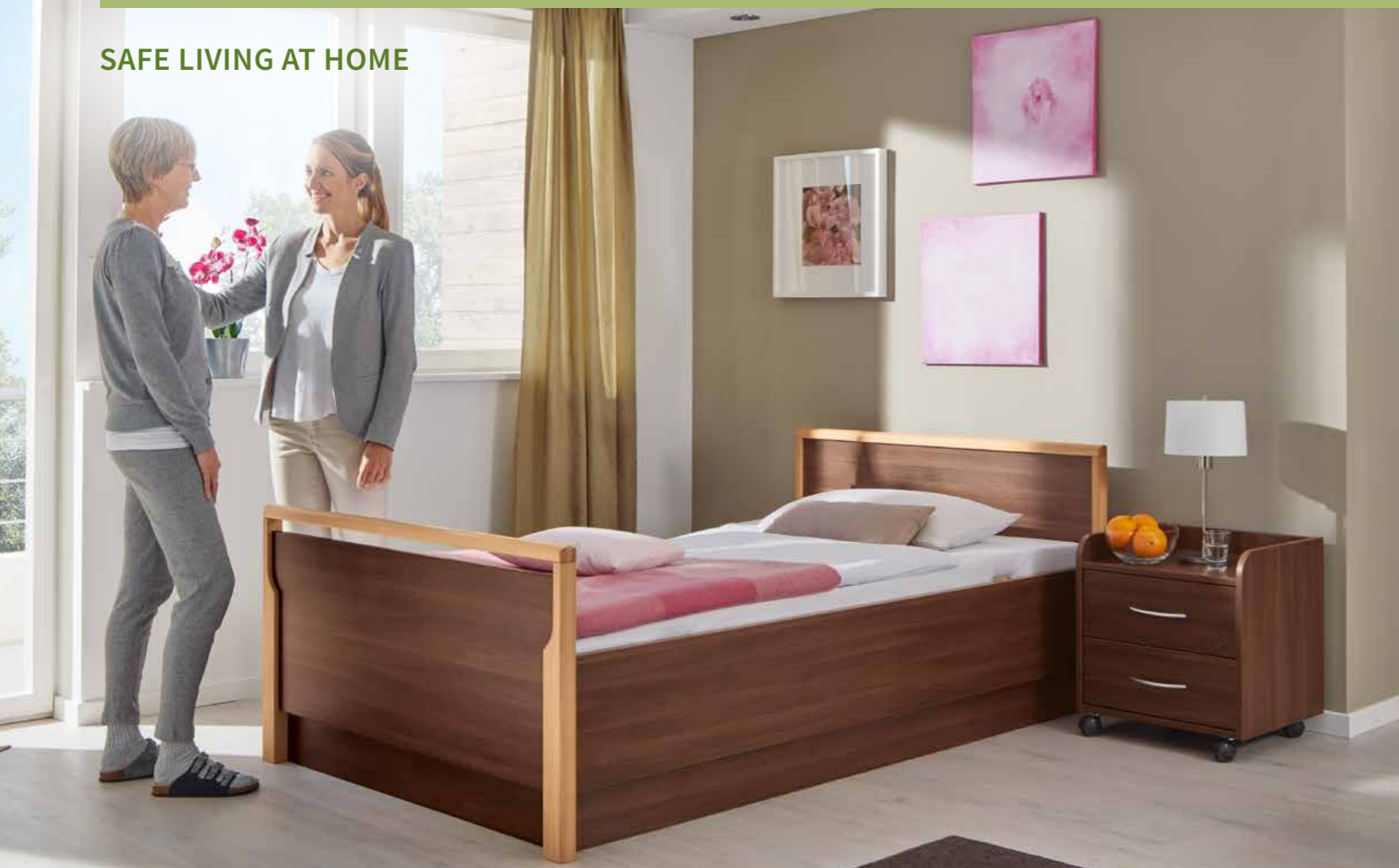

The combinable, home-comfort bed surround

Relax

ADVANTAGES

The Relax is a wooden bed surround for home comfort. It can be combined perfectly with the Lippe range models to create an elegant bed. An electrically adjustable slatted frame is also available. The Relax is also suitable for people who are looking for a double bed with technical support. Two bed surrounds can be combined to create a double bed. Each person can then choose a mattress base to suit their individual needs.

- Mattress base dimensions of 90, 100, 120 and 140 cm
- Lengths of 200 and 220 cm
- Headboard and footboard frames made of solid beech
- Panels and bed sides in Natural Beech, Royal Maple, Havana Cherry or Bella Noce Choco
- Optional wood cover to bridge the middle gap of a double bed



Further information:



Relax double bed with height adjustable Lippe IV bed insert.

The advantages at a glance

DESIGN DETAILS



Wooden frame in solid beech; panels and bed sides in Natural Beech decor.



Wooden frame in solid beech; panels and bed sides in Royal Maple decor.



Wooden frame in solid beech; panels and bed sides in Havana Cherry decor.



Wooden frame in solid beech; panels and bed sides in Bella Noce Choco decor.



Height-adjustable strips for supporting slatted frame (4 heights in 2-cm stages). Lowest position 26 cm.

ACCESSORIES



Mobilisation aid combined with Lippe IV (similar to Lippe 120 illustration).



Clamp-on Lippe safety side combined with Lippe IV.



Adjustable safety sides combined with Lippe IV.



Patient lifting pole with triangular handle combined with Lippe IV.



Wooden cover to join the gap between beds when used as a double bed.



Motorised frame; electrically adjustable backrest and leg rest (90 x 200 cm and 100 x 200 cm).



Lippe IV height-adjustable bed insert; electrically adjustable backrest, leg rest and mattress base height.

EQUIPMENT AND OPTIONS

Equipment/Option	Relax	Material no.
Design features		
Relax 90 x 200 cm		214451LFM007
Relax 100 x 200 cm		214451LFM011
Relax 120 x 200 cm		214451LFM013
Relax 120 x 220 cm		214451LFM019
Relax 140 x 200 cm		214451LFM022
Relax 140 x 220 cm		214451LFM023
Technical data		
Overall dimensions	100 x 209 cm 110 x 209 cm 130 x 209 cm 130 x 229 cm 150 x 209 cm 150 x 229 cm	
Standard features		
Mattress base size	90 x 200 cm 100 x 200 cm 120 x 200 cm 120 x 220 cm 140 x 200 cm 140 x 220 cm	
Special equipment and accessories		
Floor-level sides	O	HOL067
Slatted frame, electrically adjustable, 90 x 200 cm	O	215254
Slatted frame, electrically adjustable, 100 x 200 cm	O	215255
Adjustable safety sides for use with Lippe IV	Right Left	250138 250139
Bella Noce Choco	O	
Natural Beech	O	DEM027
Royal Maple	O	DEM022
Havana Cherry	O	DEM024
Mattress piece for bed extension, for use with Lippe 120	O	241942
Mattress piece for bed extension, for use with Lippe 140	O	244859
Reading lamp for easy insertion, for use with Lippe IV, Lippe 120/140	O	264877
Reading lamp for fixing to the patient lifting pole, for use with Lippe IV, Lippe 120/140	O	1 x 191879 1 x 174231
Transmitter holder for handset, for use with Lippe IV	O	199137
Handset extension cable, for use with Lippe IV	O	178986

Regia/Inovia



Regia/Inovia





A bed to satisfy discerning tastes


Regia

ADVANTAGES

The Regia is an elegant care bed for discerning private customers. Being a low-height bed with a large height adjustment range and optional split safety sides with the flexible Easy Switch system, it protects Regia occupants from fall injuries without the use of liberty-depriving measures.

- Safe working load: 225 kg
- Maximum weight of resident: 195 kg
- Height adjustment range: 26 to 80 cm
- 24-volt drive system
- No leakage current measurements during expected service life
- Reverse-Trendelenburg position
- Sliding backrest in accordance with the recommendations of the German Nurses Association (DBfK)
- Split telescopic safety sides with Easy Switch, or full-length safety sides
- A selection of attractive wood decors and fabric or leather upholstery

External dimensions/height adjustment range:

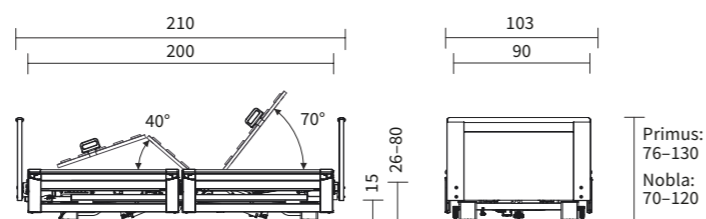

Special features:



Easy Switch

Sliding backrest in accordance with the recommendations of the German Nurses Association (DBfK)

Safety sides with 2-stage vertical telescopic extension

Electrical stepless adjustment of lower leg rest (auto-contour function)

A selection of homelike wooden surrounds

Safety sides, headboards, footboards and panels with Easy Switch system, requiring no tools for fitting

Castors with elegant covers and central locking

Optional comfort mattress base

Further information:



Product page

Assembly video

The advantages at a glance

HEAD AND FOOTBOARDS, DECORS AND UPHOLSTERY



Primus, closed surface headboard, footboard with handle bar (wood decor only).
Primus headboard, Primus-Low at foot end (wood decor only).
Nobla (upholstered only).
Nobla headboard, Nobla reha footboard (upholstered only).
Prinzino with handle bar (wood decor only).



Natural Beech (standard).
Havana Cherry.
Lindberg Oak.
Royal Maple.
Real Cherry (for Primus and Primus-Low head/footboards only).
Version with click-on full-length safety sides.

ACCESSORIES



The approximately 20-cm bed extension is simply inserted into the bed frame. Suitable mattress pieces are also available.



The storage tray for glasses, mobiles or care utensils can be attached to the full-length safety side.



The document holder is suitable for medical documents or magazines.



The comfort mattress base comprises 50 free-floating spring elements.



Wall deflection rollers protect the bed and room furnishings.



Wall spacers protect against collisions.



The under bed light ensures better orientation at night.



The linen holder is particularly practical when making beds.



Foam leather covers protect from knocks against the safety side.



Flexible transmitter holder for the handset.



Extension cable for the handset (1.40 m).



The reading lamp can be fitted to the patient lifting pole or into the connection sleeve.

EQUIPMENT AND OPTIONS

Equipment/Option	Regia	Material no.
Design features		
Regia care bed 90 x 200 cm	4-section mattress base with metal slats, chassis, undercarriage, height adjustment 26 – 80 cm, centrally locking castors, and castor covers for a homelike touch	ZLFB000
Regia care bed 100 x 200 cm	The extra cost includes the 100-cm-wide mattress base	ZLFB002
Side elements		
Split safety sides, on both sides, telescopic (TSG)	2 x TSG 110 cm, 2 x TSG 90 cm, attachment at head or foot end as required	SEE_001
Split safety sides at head end, telescopic (TSG), and wooden side panels	2 x TSG at head end 110 cm, 2 x wooden side panels at foot end 90 cm. TSG 110 cm is a fully standard-compliant safety side.	SEE_004 SEA_001 SBU_001
Split safety side at head end, telescopic (TSG), and wooden side panels	2 x TSG at head end 90 cm, 2 x wooden side panels at foot end 110 cm. TSG 90 leaves plenty room for getting in and out of bed.	SEE_004 SEA_002 SBU_001
Split safety sides at head end, telescopic (TSG), and fabric/faux leather side panels	2 x TSG at head end 110 cm, 2 x fabric/faux leather side panels at foot end 90 cm. TSG 110 cm is a fully standard-compliant safety side.	SEE_004 SEA_001 SBU_002
Split safety side at head end, telescopic (TSG), and fabric/faux leather side panels	2 x TSG at head end 90 cm, 2 x fabric/faux leather side panels at foot end 110 cm. TSG 90 leaves plenty room for getting in and out of bed.	SEE_004 SEA_002 SBU_002
Wooden side panels	2 x wooden side panels 110 cm, 2 x wooden side panels 90 cm, attachment at head or foot end as required	SEE_003 SBU_001
Fabric/faux leather side panels	2 x fabric/faux leather side panels 110 cm, 2 x fabric/faux leather side panels 90 cm, attachment at head or foot end as required	SEE_003 SBU_002
Full-length safety sides on both sides (DSG)	4 x safety side bars 200 cm, 2 x side panels 200 cm	SEE_002

Headboards and footboards		
Primus	1 x closed surface headboard with solid wood bar, 1 x footboard with open handle and solid wood bar	KFH_022
Primus-Low	1 x closed surface headboard with solid wood bar, 1 x footboard with low-level open handle and solid wood bar for a clear view of the room	KFH_023
Prinzino	1 x headboard with open handle bar made of solid wood 1 x footboard with open handle bar made of solid wood	KFH_004
Nobla (upholstered), cannot be combined with full-length safety sides (DSG)	1 x upholstered headboard, 1 x upholstered footboard. Fixed upholstery with non-removable fabric/faux leather.	KFH_019
Nobla reha (upholstered), cannot be combined with full-length safety sides (DSG)	1 x upholstered headboard, 1 x upholstered low-height footboard. Fixed upholstery with non-removable fabric/faux leather.	KFH_020
Cover materials, faux leather (The price is included when ordering upholstered headboards and side panels.)	Light green Tiba 6210 Orange Elvin 3550 Red Tiba 3210 White faux leather 1-1142-197 Anthracite faux leather 1-1142-395	BSH_H006 BSH_C001 BSH_H001 BSH_K009 BSH_K025


Technical data		
Net weight	135 kg	
Overall dimensions	110 x 210 cm, 100 x 210 cm	
Height adjustment range	26-80 cm	
Safe working load	225 kg	
Maximum resident weight	195 kg	
24-volt system	✓	

Standard features		
Frame colour	Argentum	
Decor	Natural Beech	
Mattress base size	140 x 200 cm, 120 x 200 cm	
Metal slatted mattress base (mattress base model)	✓	
Ratchet mechanism, electrical adjustment of lower leg rest for metal slatted mattress base	✓	

Select fixed upholstery with attractive fabrics and faux leather for the head/footboards and side panels of the Regia.

Fabrics:
Light green Tiba 6210
Orange Elvin 3550
Red Tiba 3210

Faux leather:
White 1-1142-197
Anthracite 1-1142-395

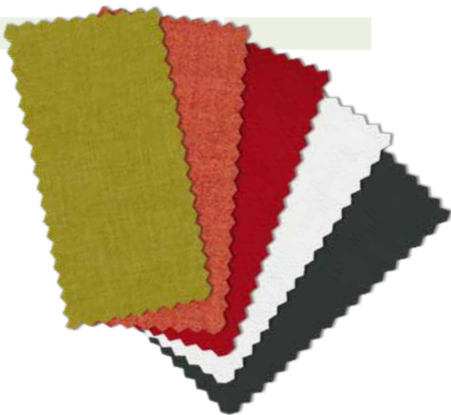

✓ = Standard feature O = Available as an option

Equipment/Option	Regia	Material no.
Special equipment and accessories		
Patient lifting pole	O	AUI_001
Havana Cherry	O	DHU_024
Real Cherry (for Primus and Primus-Low head/footboards only)	O	DHU_012
Royal Maple	O	DHU_022
Lindberg Oak	O	DHU_044
Bed extension for metal slatted mattress base 90 and 100 cm with TSG	O	LFL_008
Bed extension for metal slatted mattress base with DSG	O	LFL_005
Mattress piece for bed extension 90 cm wide, 14 cm thick	O	209240
Mattress piece for bed extension 100 cm wide, 14 cm thick	O	250270
Linen holder	O	IBA_001
Trendelenburg handset	O	ZHS_003
Comfort mattress base	O	ZLFA003
Wall spacer, 1 piece, head end/side of bed	O	262024/ 262000
Foam leather cover for full-length safety sides (DSG)	O	225921
Foam leather cover for telescopic safety sides (TSG)	O	716556
Reading lamp for easy insertion	O	264877
Reading lamp for fixing to the patient lifting pole	O	1 x 191879 1 x 174231
Under bed light	O	UBB_001
Tray for full-length safety sides (DSG), 90 cm wide	O	174296
Tray for full-length safety sides (DSG), 100 cm wide	O	273465
Extension cushion (for full-length safety side only)	O	258395
Transmitter holder for handset	O	199137
Handset extension cable	O	238047
Document holder	O	204241
Carry bag for safety side bars (DSG)	O	273831
Storage tray	White Yellow Red	266970 266971 266619



The care bed with innovative ideas


Inovia

ADVANTAGES


The Inovia care bed has an attractive, homelike appearance thanks, in particular, to its opulent headboards and footboards. Its lowest position, at just 22 cm, protects resident from falls, without having to resort to liberty-depriving measures. The optional telescopic safety sides also promote the resident's mobility.

- Safe working load: 225 kg
- Maximum weight of resident: 195 kg
- Height adjustment range: 22 to 77 cm
- 24-volt drive system
- No leakage current measurements during expected service life
- Reverse-Trendelenburg position
- Sliding backrest in accordance with the recommendations of the German Nurses Association (DBfK)
- Optional telescopic split safety side
- Attractive wood decors to choose from

External dimensions/height adjustment range:

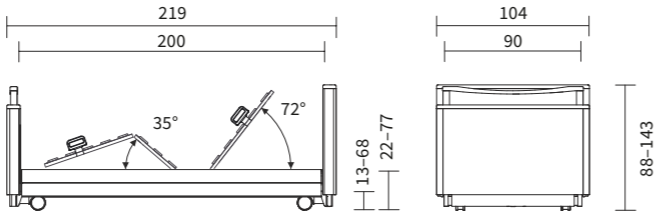

Special features:

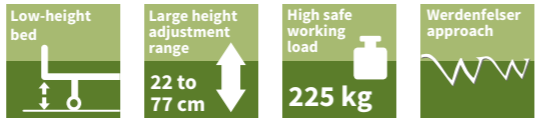

Further information:



The advantages at a glance

DESIGN DETAILS



Solid wood handle bar.



Castors can be locked in pairs whatever the mattress base position.



Handset with selective locking.



Sliding backrest.



Mattress base height adjustable from 22 to 77 cm, adjustment range = 55 cm.



Natural Beech (standard).



Havana Cherry.

ACCESSORIES



Telescopic safety side (TSG).



Foam leather cover for DSG.



Foam leather cover for TSG.



Bed extension.



Tray.



Sola reading lamp.



The storage tray for glasses, mobiles or care utensils can be attached to the full-length safety side.



The comfort mattress base comprises 50 free-floating spring elements.



Flexible transmitter holder for the handset.

EQUIPMENT AND OPTIONS

Equipment/Option	Inovia	Material no.
Design features		
Inovia 90 x 200 cm with full-length safety sides (DSG)		262626LFB_001 SHC_001
Inovia 100 x 200 cm with full-length safety sides (DSG)		262626LFB_002 SHC_001
Technical data		
Net weight	153/159 kg	
Overall dimensions	105 x 220 cm 115 x 220 cm	
Height adjustment range	22 – 77 cm	
Safe working load	225 kg	
Maximum resident weight	195 kg	
24-volt system	✓	
Standard features		
Frame colour	Argentum	
Decor	Natural Beech	
Mattress base size	90 x 200 cm 100 x 200 cm	
Metal slatted mattress base (mattress base model)	✓	
Type of handset (selectively lockable/wireless)	Selective	
Mechanical CPR release	✓	
Castor diameter	100 mm	
Castor locking	Per axle	
Ratchet mechanism, electrical adjustment of lower leg rest	✓	
Mechanical CPR release	✓	
Special equipment and accessories		
Telescopic safety sides (TSG) on both sides	O	SHC_002
Telescopic safety side (TSG) on one side/DSG on one side	O	SHC_008
Patient lifting pole	O	AUI_001
Natural Beech	✓	DHU_027
Havana Cherry	O	DHU_024
Handle bar at head end	O	KTA_003
Bed extension for metal slatted mattress base 90 and 100 cm with TSG	O	262626LFL_006
Bed extension for metal slatted mattress base 90 and 100 cm with DSG	O	262626LFL_005
Mattress piece for bed extension 90 cm wide, 12 cm thick	O	166715
Mattress piece for bed extension 100 cm wide, 12 cm thick	O	219857
Comfort mattress base	O	LIF_013
Wall deflection roller with holder (right or left)	O	Left 215921 Right 221306
Foam leather cover for full-length safety sides (DSG)	O	225921
Foam leather cover for telescopic safety sides (TSG)	O	716556
Reading lamp for easy insertion	O	264877
Reading lamp for fixing to the patient lifting pole	O	1x191879 1x174231
Tray for full-length safety sides (DSG), 90 cm wide	O	174296
Tray for full-length safety sides (DSG), 100 cm wide	O	229222
Extension cushion (for full-length safety side only)	O	258395
Transmitter holder for handset	O	199137
Handset extension cable	O	238047
Storage tray	White Yellow Red	266970 266971 266619

✓ = Standard feature O = Available as an option

Westfalia range



Westfalia range





The care bed for a more
homelike ambiance


Westfalia IV

ADVANTAGES

The Westfalia IV combines excellent comfort, advanced safety and a particularly stylish design. The height can be adjusted with the help of a stable yet sleek and elegant scissor mechanism. The headboards and footboards are available in various attractive forms.

- Safe working load: 185 kg
- Maximum weight of resident: 155 kg
- Height adjustment range: 35 to 80 cm
- 24-volt drive system
- No leakage current measurements during expected service life
- Selectively lockable handset
- Attractive head and footboard designs
- Different mattress base sizes to choose from

External dimensions/height adjustment range:

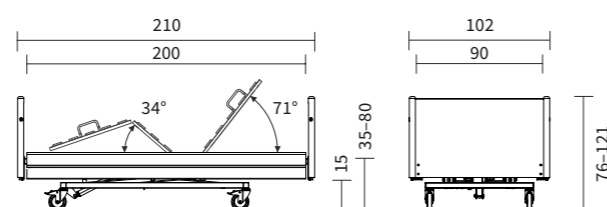

Special features:



Further information:



Product page

Assembly video

The advantages at a glance

DESIGN DETAILS



Handset with selective locking function.



Individually braked castors.



Height adjustment range 35 to 80 cm.



Patient lifting pole.



Choice of headboards and footboards: Standard (standard feature), Comodo with round posts (optional), Primero with handle bar (optional).

ACCESSORIES



Tray.



Mobilisation aid.



Adaptable safety sides.



Wall deflection roller including mount.



Reading lamps.



Bed extension including safety side bars (also for retrofitting).



Foam leather cover.



The storage tray for glasses, mobiles or the wireless handset can be attached at various positions.



Flexible transmitter holder for the handset.



Extension cable for the handset (1.40 m).



With the floor-mounted power supply unit, the 24-volt system can also be used with various international mains cables.

EQUIPMENT AND OPTIONS

Equipment/Option	Westfalia IV	Material no.
Design features		
Westfalia IV 90 x 200 cm		229911ZLFB000
Westfalia IV 100 x 200 cm ¹⁾		229911ZLFB002
Westfalia IV 90 x 180 cm		229911ZLFL001
Westfalia IV 90 x 190 cm		229911ZLFL010
Technical data		
Net weight	136 kg	
Overall dimensions	110 x 232	
	110 x 212	
	100 x 232	
	100 x 192	
	100 x 202	
	100 x 212	
	90 x 212	
Height adjustment range	35 – 80 cm	
Safe working load	185 kg	
Maximum resident weight	155 kg	
24-volt system	✓	
Standard features		
Frame colour	Topaz	
Decor	Natural Beech	
Mattress base size	90 x 200 cm	
Sprung wooden slat mattress base (mattress base model)	✓	
Safety side	Full-length safety side (DSG)	
Type of handset (selectively lockable/wireless)	Selective	
Mechanical CPR release	✓	
Castor diameter	125 mm	
Castor locking	Single	
Ratchet mechanism, electrical adjustment of lower leg rest	✓	
Patient lifting pole	✓	
Special equipment and accessories		
Adaptable safety sides	Right	250138
	Left	250139
Mobilisation aid	O	175402
Comodo wooden surround with round posts	O	ZMOK015
Primero wooden surround	O	ZMOK017
Bed extension to mattress base length of 220 cm	O	ZLFL007
Mattress piece for bed extension 90 cm wide, 12 cm thick	O	166715
Mattress piece for bed extension 100 cm wide, 12 cm thick	O	219857
4 wall deflection rollers, horizontal action	O	ZWAN004
2 wall deflection rollers, horizontal action, foot end	O	ZWAN003
2 wall deflection rollers, horizontal action, head end	O	ZWAN002
Foam leather cover for full-length safety sides (DSG)	O	225921
Reading lamp for easy insertion	O	264877
Reading lamp for fixing to the patient lifting pole	O	1 x 191879
		1 x 174231
Tray for full-length safety sides (DSG), 90 cm wide	O	174296
Tray for full-length safety sides (DSG), 100 cm wide	O	273465
Extension cushion (for full-length safety side only)	O	258395
Transmitter holder for handset	O	199137
Handset extension cable	O	178986
Document holder	O	204241
Storage tray	White	266970
	Yellow	266971
	Red	266619
Floor-mounted 24-volt power supply unit	O	254276

✓ = Standard feature O = Available as an option ¹⁾ With Comodo head/footboard only



The care bed
with special features


Westfalia-Klassik

ADVANTAGES


The Westfalia-Klassik is designed for discerning users who are not prepared to compromise on comfort, individuality and sophisticated technical features in their daily care routine.

- Safe working load: 200 kg
- Maximum weight of resident: 170 kg
- Height adjustment range: 35 to 80 cm
- 24-volt drive system
- No leakage current measurements during expected service life
- Sliding backrest in accordance with the recommendations of the German Nurses Association (DBfK)
- Central castor locking

External dimensions/height adjustment range:

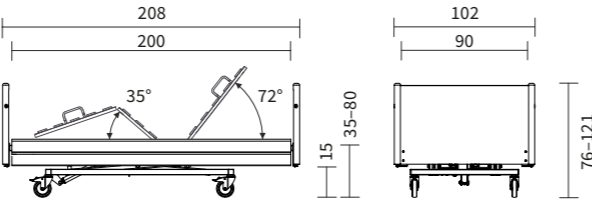

Special features:

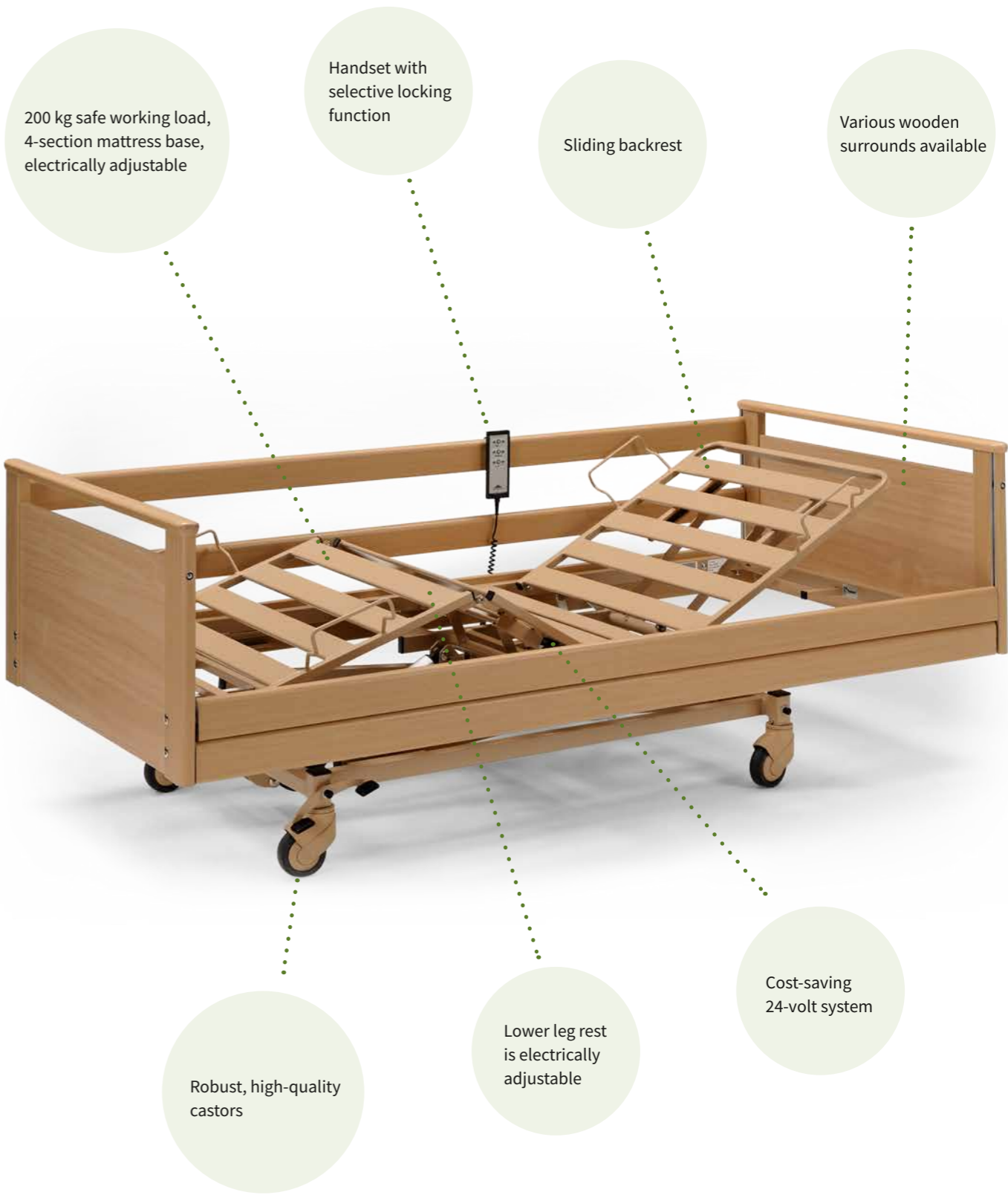

Further information:



The advantages at a glance

DESIGN DETAILS



Central braking possible whatever the mattress base position.



Electronic (SwitchMode) power supply unit integrated immediately downstream of the mains plug.



Mattress base height adjustable from 35 to 80 cm, adjustment range = 45 cm.

ACCESSORIES



Tray.




Wall deflection roller.



Wooden surrounds: Standard (standard feature), Comodo with round posts (optional), Primero with handle bar (optional).



Reading lamps.

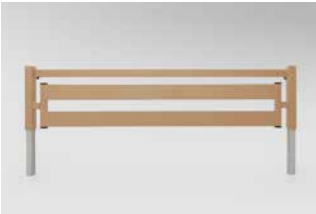

Bed extension with metal slats, safety sides and side panels.



Foam leather cover.



The storage tray for glasses, mobiles or the wireless handset can be attached at various positions.



Flexible transmitter holder for the handset.



Extension cable for the handset (1.40 m).



Extension cushion to place over the full-length safety side.



With the floor-mounted power supply unit, the 24-volt system can also be used with various international mains cables.

EQUIPMENT AND OPTIONS

Equipment/Option	Westfalia-Klassik	Material no.
Design features		
Westfalia-Klassik 90x200 cm with standard wooden surround		171358HOL006
Technical data		
Net weight	115 kg	
External dimensions with standard wooden surround	102 x 208 cm	
Height adjustment range	35 – 80 cm	
Safe working load	200 kg	
Maximum resident weight	170 kg	
24-volt system	✓	
Standard features		
Frame colour	Topaz	
Decor	Natural Beech	
Mattress base size	90 x 200 cm	
Safety side	Full-length safety side (DSG)	
Type of handset (selectively lockable / wireless)	Selective	
Mechanical CPR release	✓	
Castor diameter	100 mm	
Castor locking	Central	
Ratchet mechanism, electrical adjustment of lower leg rest	✓	
Accessories		
Comodo wooden surround with round posts	O	HOL007
Primero wooden surround	O	HOL012
Bed extension	O	701215 + 203806
Wall deflection roller Ø 100 mm	O	WAN011
Foam leather cover for full-length safety sides (DSG)	O	225921
Reading lamp for easy insertion	O	264877
Reading lamp for fixing to the patient lifting pole	O	1 x 191879 1 x 174231
Tray for full-length safety sides (DSG)	O	174296
Extension cushion (for full-length safety side only)	O	258395
Transmitter holder for handset	O	199137
Handset extension cable	O	178986
Document holder	O	204241
Storage tray	White	266970
	Yellow	266971
	Red	266619
Floor-mounted 24-volt power supply unit	O	204276

✓ = Standard feature O = Available as an option

Heavy-duty beds



Heavy-duty beds





The heavy-duty care bed


Gigant

ADVANTAGES

The Gigant heavy-duty bed can support occupants weighing up to 280 kg, making it suitable for extremely obese persons. It also offers plenty of space and support for independent living and reduces the burden on carers.

- Safe working load: 350 kg
- Maximum weight of resident: 280 kg
- Height adjustment range: 40 to 80 cm
- Reverse-Trendelenburg position
- Optional combinable safety side
- Headboards and footboard with elegant wooden surround

External dimensions/height adjustment range:

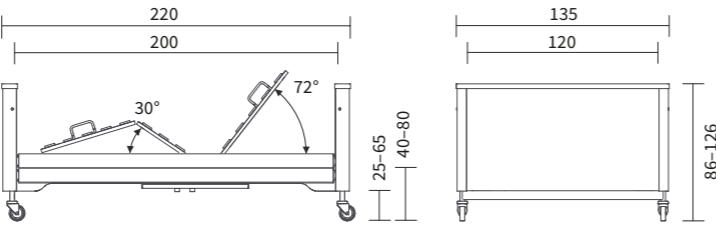

Special features:



Further information:



DESIGN DETAILS



Holder for mobilisation aid (left) and locking box (right).



CPR release of the backrest.



Extremely sturdy safety side construction.



Plastic panels on the chassis for a more homelike feel.

ACCESSORIES



Clamp-on safety side bar.



Combinable safety side.



Reading lamps.



Bed extension including longer safety sides and side panels.



Foam leather cover for telescopic safety sides (TSG).




Extension cushion to place over the full-length safety side.



The storage tray for glasses, mobiles or the wireless handset can be attached at various positions.

EQUIPMENT AND OPTIONS

Equipment/Option	Gigant	Material no.
Design features		
Gigant 120 x 200 cm		238647LFM013
Technical data		
Net weight	159 kg	
Overall dimensions	135 x 222 cm	
Height adjustment range	40 – 80 cm	
Safe working load	350 kg	
Maximum resident weight	280 kg	
Standard features		
Frame colour	Argentum	
Decor	Natural Beech	
Mattress base size	120 x 200 cm	
Safety side	Full-length safety side (DSG)	
Type of handset (selectively lockable / wireless)	Selective	
Mechanical CPR release	✓	
Castor diameter	100 mm	
Castor locking	Single	
Ratchet mechanism, electrical adjustment of lower leg rest	✓	
Accessories		
Clamp-on safety side bar	O	164977
Combinable safety side on both sides	O	SEG022
Combinable safety side on left side	O	SEG023
Combinable safety side on right side	O	SEG024
Bed extension for metal slatted mattress base (with DSG)	O	240271 + 242308
Mattress piece for bed extension 120 cm wide, 19 cm thick	O	712167
Foam leather cover for DSG, one side	O	225921
Foam leather cover for KSG, one side	O	2x 183008 + 1x 712961
Reading lamp for easy insertion	O	264877
Reading lamp for fixing to the patient lifting pole	O	1x 191879 1x 174231
Extension cushion (for full-length safety side only)	O	258395
Storage tray	White Yellow Red	266970 266971 266619

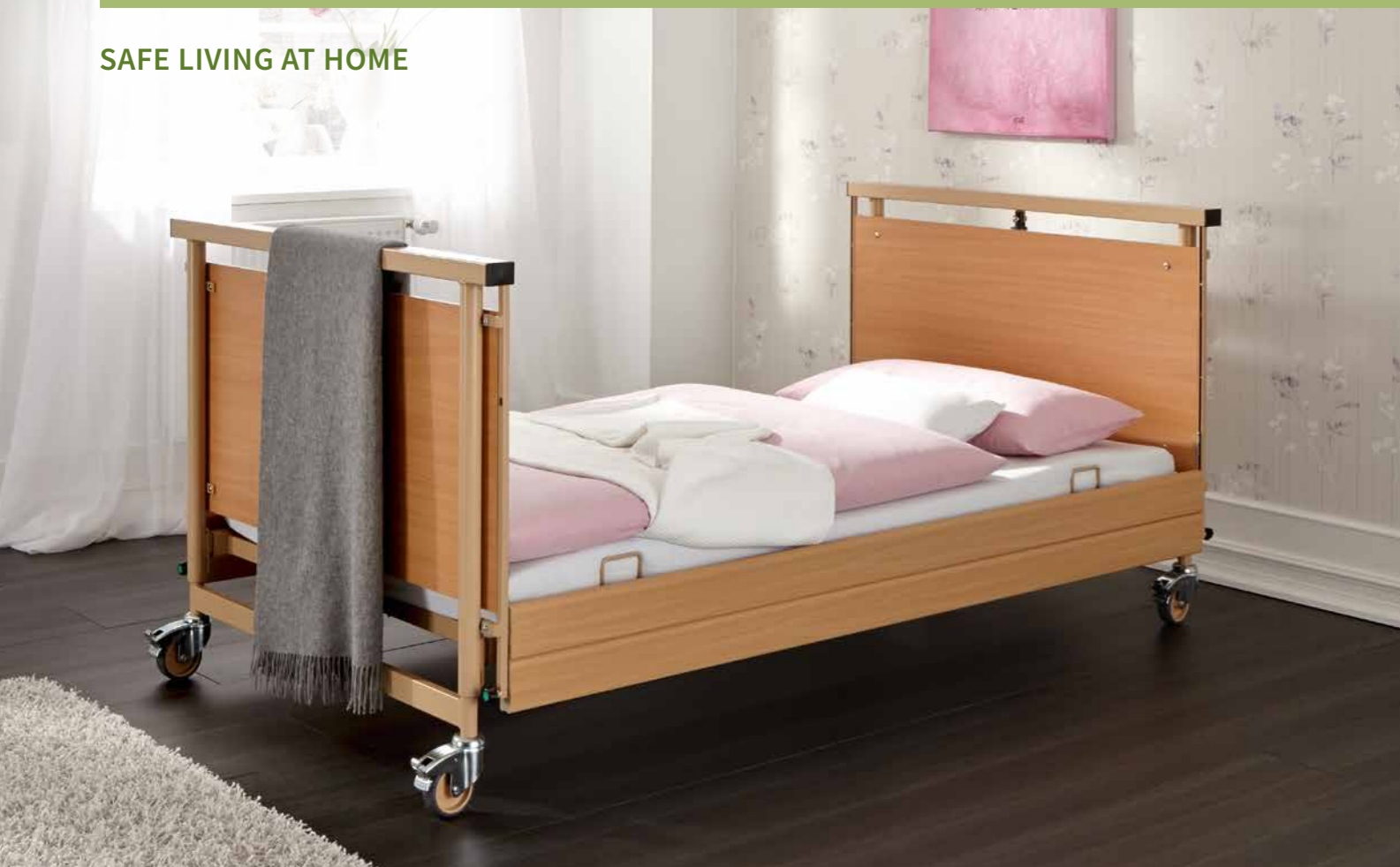

The innovative
care beds


Allura II

ADVANTAGES

The flexible models of the Allura II range all have a high load capacity with a choice of mattress base widths and a large height adjustment range. This means that they offer even heavy residents a high level of comfort and safety, while relieving the burden placed on carers.

- Safe working load: 250 kg
- Maximum weight of resident: 210 kg
- Height adjustment range: 30 to 80 cm
- 24-volt drive system
- No leakage current measurements during expected service life
- Reverse-Trendelenburg position
- Selectively lockable handset
- Mattress base widths of 100 and 120 cm

External dimensions/height adjustment range:

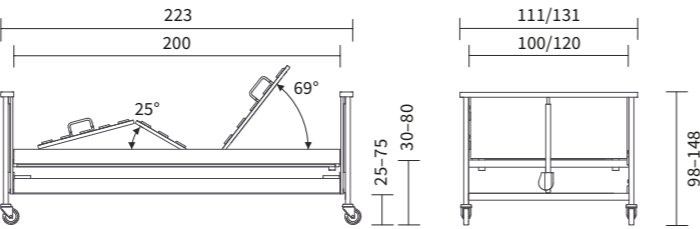


Special features:

Heavy-duty bed

High safe working load
250 kg

German engineering design

24 volt
system

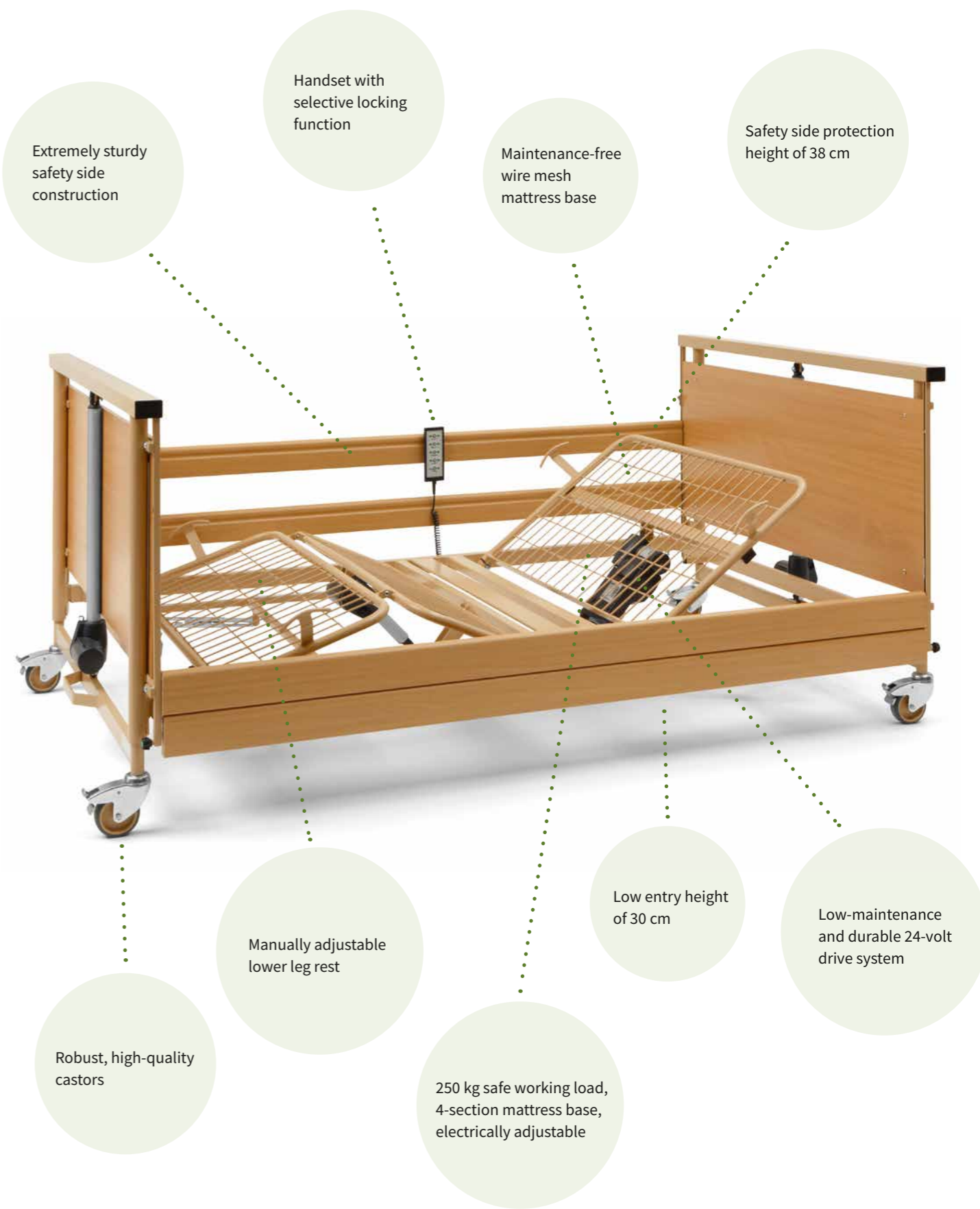

Further information:



The advantages at a glance

DESIGN DETAILS



Reverse-Trendelenburg position.



Adjustment range from 30 to 80 cm.



Electronic (SwitchMode) power supply unit integrated immediately downstream of the mains plug.



Handset with selective locking function.



Compact delivery and storage on storage aid/transportation aid (similar to illustration).

ACCESSORIES



Bed extension (also for retrofitting).



Foam leather cover.



Reading lamps.



Tray for Allura II 100.



Wooden panels to cover the chassis (also for retrofitting).



Extension cushion to place over the full-length safety side.



The storage tray for glasses, mobiles or the wireless handset can be attached at various positions.



Wall deflection roller including mount.



The document holder is suitable for medical documents or magazines.

EQUIPMENT AND OPTIONS

Equipment/Option	Allura II	Material no.
Design features		
Allura II 100x200 cm		275100
Allura II 120x200 cm		275101
Allura 100x220 cm		278115
Allura 120x220 cm		278116
Technical data/standard features		
Net weight	120 kg	
Overall dimensions	111 x 224 cm	
Height adjustment range	30 – 80 cm	
Safe working load	250 kg	
Maximum resident weight	210 kg	
24-volt system	✓	
Standard features		
Frame colour	Topaz	
Decor	Natural Beech	
Mattress base size	100 x 200 cm	
Safety side	Full-length safety side (DSG)	
Type of handset (selectively lockable/wireless)	Selective	
Mechanical CPR release	✓	
Castor diameter	100 mm	
Castor locking	Single	
Ratchet mechanism, electrical adjustment of lower leg rest	✓	
Accessories		
Mobilisation aid	O	244029
Bed extension for metal slatted mattress base 100 cm	O	1x 195915 4x 198167
Bed extension for metal slatted mattress base 120 cm	O	1x 195906 4x 198167
Mattress piece for bed extension 100 cm wide, 12 cm thick	O	219857
Mattress piece for bed extension 120 cm wide, 16 cm thick	O	241942
Wall deflection roller, 1 piece, including mount	O	206751
Foam leather cover for full-length safety sides (DSG)	O	225921
Reading lamp for easy insertion	O	264877
Reading lamp for fixing to the patient lifting pole	O	1 x 191879 1 x 174231
Tray for full-length safety sides (DSG) (for 100 cm width only)	O	180039
Extension cushion (for full-length safety side only)	O	258395
Document holder	O	204241
Storage tray	White Yellow Red	266970 266971 266619

✓ = Standard feature O = Available as an option

Vertica homedcare



Vertica homecare





Gentle
mobilisation


Vertica homecare

ADVANTAGES

At the touch of a button, the Vertica homecare mobilisation bed moves the resident into a genuine sitting position or an ergonomic stand-up position. It follows the body's natural sequence of movements. People in need of care are able to participate in conversations at eye level again and enjoy mealtimes. The bed helps residents to get walking again after a fall or operation.

- Safe working load: 225 kg
- Maximum weight of resident: 195 kg
- Height adjustment range: 40 to 80 cm
- Sliding backrest in accordance with the recommendations of the German Nurses Association (DBfK)
- Separate carers' handset for the sitting and stand-up positions
- 2/3 safety sides, MultiFlex
- Lindberg Oak homelike decor

External dimensions/height adjustment range:

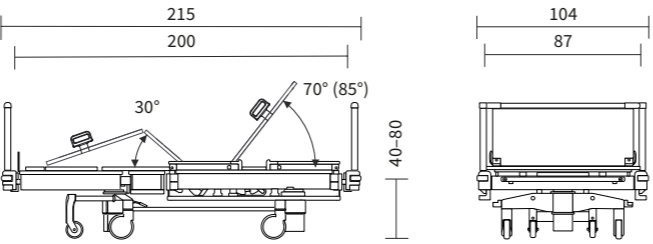

Special features:

Mobilisation 	Large height adjustment range 40 to 80 cm 	High safe working load 225 kg 	Trendelenburg position 	Manoeuvrability
-------------------------	---	---	-----------------------------------	----------------------------



Further information:



- Product video Vertica video

The advantages at a glance

DESIGN DETAILS



Separate handset for the sitting and stand-up positions.



Main handset with easily understandable pictograms.



Stable mobilisation supports.



Integrated handles provide additional support.

ACCESSORIES



The fold-down footboard can also be used as a table.



The footboard can also be used as a linen holder.



Filler piece used for extending the bed.



Mattress piece for extending the mattress.



Perfectly fitting viscoelastic mattress.



Footrests promote spatial perception.



Side support cushions give additional stability.

EQUIPMENT AND OPTIONS

Equipment/Option	Vertica homecare	Material no.
Design features		
Vertica homecare with Lindberg Oak decor		247 348
Technical data		
Net weight	Approx. 175 kg	
Overall dimensions	103.8 x 214.2 cm	
Height adjustment range	40 – 80 cm	
Safe working load	225 kg	
Maximum resident weight	195 kg	
Standard features		
Frame colour	Argentum	
Decor	Lindberg Oak R4223 W	
Mattress base size	87 x 200 cm	
Fully motorised drive	✓	
Battery (not with sit-up and stand-up function)	✓	
Polypropylene hygiene mattress base for belt system	✓	
Bed extension 15 cm	✓	
Multiflex ¾ safety sides	✓	
4 wall deflection rollers, all sides	✓	
Integral castor, 125 mm	✓	
Electrically conductive castor: 1 castor (directional wheel)	✓	
Brake pedal: cast iron	✓	
Handset with locking function	✓	
Vertica handset	✓	
Brake signal	✓	
Parking position for footrest (14 cm)	✓	
Control unit cover	✓	
Integrated bed extension	✓	
Special equipment and accessories		
Fold-down footboard for use as table and linen holder	O	240 245
Footrests	O	240 434
Special mattress with cut-out sections	O	247948
Side support cushion as armrest, left	O	709575
Side support cushion as armrest, right	O	709578
Extension piece for extending the mattress base +	O	706179 +
Mattress piece for extending the mattress (including extension board)	O	209240

✓ = Standard feature O = Available as an option

Bedside cabinets and side tables



Bedside cabinets and side tables





Practical and functional

Bedside cabinets

Further information:



Bedside cabinets and side tables



Item number: 193753

Cherusker bedside cabinet

With height-adjustable tilting table, plastic-coated body, beech decor

Technical data:

- Width: 52 cm
- Height: 87 cm
- Overall depth: 50 cm
- Overall width (with table horizontal): 114 cm
- Table height adjustment from 82 to 112 cm
- Double castors Ø 50 mm, 2 castors are lockable



Item number: 236959

Arminius bedside cabinet

With height-adjustable tilting table, plastic-coated body, beech decor

Technical data:

- Width: 54.5 cm
- Height: 89 cm
- Overall depth: 45 cm
- Overall width (with table horizontal): 125 cm
- Table height adjustment from 78 to 112 cm
- Double castors Ø 50 mm, 2 castors are lockable





Item number: 223636

Hermann bedside cabinet with table

With fold-down table, plastic-coated body, beech decor

Technical data:

Width: 55 cm
Height: 89 cm
Overall depth: 45 cm
Overall width (with table horizontal): 119.5 cm
Double castors Ø 50 mm, 2 castors are lockable



Item number: 243335

Hermann bedside cabinet without table

Plastic-coated body, beech decor, no table

Technical data:

Width: 44.5 cm
Height: 89 cm
Overall depth: 45 cm
Double castors Ø 50 mm, 2 castors are lockable



Item number: 232644

Varus bedside cabinet

Plastic-coated body, no table, available in four different decors: Natural Beech, Royal Maple, Cherry Havana, Bella Noce Choco

Technical data:

Width: 54 cm
Height: 57 cm
Overall depth: 45 cm
Double castors Ø 50 mm, 2 drawers



Item number: 218931

Duo overbed table and side table

With height-adjustable, tilting overbed table and fixed table section, beech decor

Technical data:

Width: 81 cm
Height: 81 cm
Overall depth: 42 cm
Overbed table height adjustable from 80 to 113 cm
Double castors Ø 50 mm, 2 castors are lockable



Item number: 218932

Solo overbed table and side table

With height-adjustable tilting table, beech decor

Technical data:

Width: 60 cm
Height: 81 cm
Overall depth: 42 cm
Overbed table height adjustable from 80 to 113 cm
Double castors Ø 50 mm, 2 castors are lockable

*Transport chair/
mattresses/furniture*



Transport chair, mattresses, furniture





The care chair for extra mobility

Ravello

ADVANTAGES

The Ravello transport chair allows persons in need of care to be transferred safely from the bedroom to the bathroom or dining room, for example. Its flexible construction makes it easy for carers to transfer the resident to the chair from an armchair or the side of the bed. The cushioned, ergonomically shaped seat also makes the Ravello comfortable to sit in for longer periods. The washable seat covers are very easy to keep clean.

- Safe working load of 135 kg (dynamic), 200 kg (static)
- The lowerable armrests make it easy to get in to the chair from the side.
- The backrest can be lowered to a 43° reclining position.
- Ravello-sano: retractable footrest
- Ravello-curo: permanently fixed footrest moves with the lower leg rest during adjustments
- Ravello sano-s: 4 steerable castors for easy manoeuvring
- Ravello plus-s: height-adjustable seat

Further information:



Product page



The advantages at a glance

DESIGN DETAILS



Ravello sano: the footrest can be completely retracted.



Ravello curo: length-adjustable, fixed position footrest that tilts to follow backrest adjustments.



Armrest height adjustment in 6 stages.



The Ravello sano-s is available for better manoeuvrability.

ACCESSORIES



Side supports.



Headrest with lateral support.



Urine bag holder.



Attachable table.



Infusion stand.

EQUIPMENT AND OPTIONS

Equipment/Option	Ravello	Material no.
Design features		
Ravello curo		264299
Ravello sano		264297
Ravello sano-s		264298
Ravello plus-s		264300
Technical data		
Net weight	Curo, Sano, Sano-s: 38 kg Plus-s: 58 kg	
Dimensions	Curo, Sano, Sano-s: W x H x D 68 x 126 x 103 cm incl. footrests Plus-s: W x H x D 73 x 132–161 x 109 cm incl. footrests	
Backrest	Height 81 cm, reclines to 43°	
Footrest	26 x 24 cm, with 8 cm ground clearance	
Castors	Front castors: 20 cm, fixed Rear castors: 15 cm, swivel type, with paired locking	
S version castors:	15 cm, swivel type, with paired locking	
Seat	PU foam with wooden core	
Backrest	PU foam with steel core	
Covering fabric	Faux leather	
Special equipment and accessories		
Tray, attachable table	O	195900
Neck support (black)	O	116344
1 pair of side supports (azure blue)	O	189619
Leg support, right	O	116342
Leg support, left	O	116343
Urine bag holder	O	116340
Infusion stand	O	1x 260269 1x 114352

✓ = Standard feature O = Available as an option



Well-being through comfort

Mattresses

ADVANTAGES


Enjoy a good night’s sleep
Burmeier’s care beds have a range of adjustment options for comfort and support, day and night. Just as important for restful sleep, of course, is the choice of the right mattress. Our mattresses are ideally suited to the various requirements of care in the home. They ensure a healthy, dry microclimate, adapt flexibly to the shape of the mattress base and reduce pressure point loading. The bed and mattress can be selected as a perfectly matching unit.

Pressure ulcer protection
People in need of care must be effectively protected from pressure point loading during prolonged periods in bed. All Burmeier mattresses make an important contribution here. They continuously adapt to the shape of the body and ensure uniform pressure distribution.

Further information:



Mattresses for a new sleep sensation

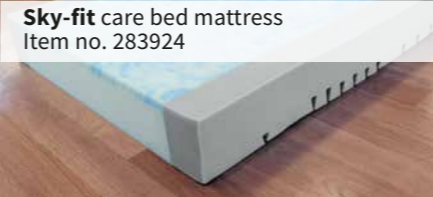

The premium model to satisfy every requirement

- The high-quality dual-layer mattress for Burmeier beds has a 2 cm Sky-therm® layer and a 12 cm base layer made of cold foam
- Reinforced edge zones
- Adapts perfectly to the mattress base and follows every bed adjustment
- Removable, washable and reversible cover; one side is PU and the other polyester and viscose
- Compression hardness 5.0 kPa (edge zones), 4.0 kPa (base layer), 2.5 kPa (top layer), density 50 kg/m³ (base), 40 kg/m³ (top)



The high-quality mid-range mattress for healthy sleep

- The comfortable mattress with a core made of elastic polyurethane cold foam (14 cm) for recuperative sleep
- Reinforced edge zones
- Adapts perfectly to the mattress base and follows every bed adjustment
- Removable, washable and reversible cover; one side is PU and the other polyester and viscose
- Compression hardness 5.0 kPa (edge zones), 4.0 kPa (base), density 40 kg/m³



The versatile entry-level model for everyday care

- The modern care bed mattress with a polyurethane core (12 cm)
- Washable PU cover, white, removable with zip fastening
- Density: 40 kg/m³, compression hardness: 4.5 kPa

DIMENSIONS

Model	Dimensions (W x L x H)	Item no.
Sky-fit care bed mattress	90 cm x 200 cm x 14 cm	283924
	100 cm x 200 cm x 14 cm	283925
Comfort-fit care bed mattress	90 cm x 200 cm x 14 cm	283922
	100 cm x 200 cm x 14 cm	283923
Basic-fit care bed mattress	80 cm x 180 cm x 12 cm	261478
	80 cm x 190 cm x 12 cm	261479
	80 cm x 200 cm x 12 cm	261481
	90 cm x 180 cm x 12 cm	261482
	90 cm x 190 cm x 12 cm	261483
	90 cm x 200 cm x 12 cm	261477
	90 cm x 220 cm x 12 cm	261484
	100 cm x 190 cm x 12 cm	261485
	100 cm x 200 cm x 12 cm	261486
	100 cm x 220 cm x 12 cm	261487
	120 cm x 200 cm x 16 cm	261488
	140 cm x 200 cm x 16 cm	261489



Flexible and functional

Wardrobes and cabinets

ADVANTAGES

Flexible and functional
With wardrobes and cupboards from our Adrano furniture system, bedrooms and living rooms can be designed with individual style. Clear lines and a wood decor create a pleasant, cosy atmosphere. Clothes, crockery or care utensils can be stored for easy access.

Durable
The entire furniture range is manufactured with ultra-modern equipment and high quality materials. The rear panels, for example, are 8 mm thick. All furniture units have a wipe-resistant base to protect them from moisture. The underside of the base is additionally protected with ABS edging.

Further information:



Wardrobes/linen cupboards, sideboards




Adrano
offers a variety of different options for keeping valuables safe. All wardrobes and cupboards can be fitted with optional safes. These lockable compartments can also be built in as permanent fixtures.


Mix and match
Take advantage of the wide range of furniture units available to create your very own wardrobe and cabinet combination.



Examples of combinations

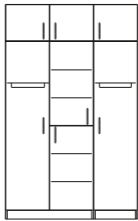

Combination of 241749 & 241755
Technical data:
Width: 150 cm
Height: 244 cm
Overall depth: 60 cm

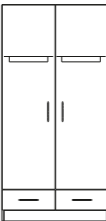


Combination of 241748 & 241739
Technical data:
Width: 150 cm
Height: 244 cm
Overall depth: 60 cm
With 3 base drawers




Material no.: 241634*
Technical data:
Width: 100 cm
Height: 90 cm
Depth: 52 cm
With 2 drawers




Material no.: 241661*
Technical data:
Width: 100 cm
Height: 90 cm
Depth: 52 cm

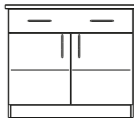

Material no.: 241662*
Technical data:
Width: 100 cm
Height: 90 cm
Depth: 52 cm
With 1 drawer

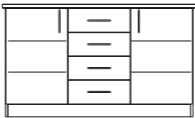

Material no.: 241633*
Technical data:
Width: 150 cm
Height: 90 cm
Depth: 52 cm
With 4 drawers

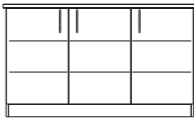


Material no.: 241664*
Technical data:
Width: 150 cm
Height: 90 cm
Depth: 52 cm


*If sideboards or cupboards are combined, the adjoining top surface of the sideboard ends flush with the sideboard body.

Wardrobes/linen cupboards, top units and accessories


One-person wardrobes/linen cupboards*

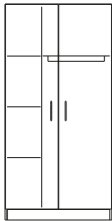


Material no.: 241740
Technical data:
Width: 100 cm
Height: 202 cm
Depth: 60 cm

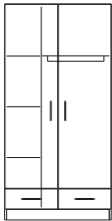

Material no.: 241741
Technical data:
Width: 100 cm
Height: 202 cm
Depth: 60 cm
With 2 base drawers

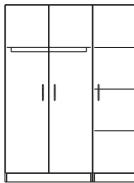


Material no.: 241742
Technical data:
Width: 150 cm
Height: 202 cm
Depth: 60 cm

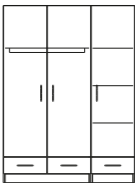


Material no.: 241743
Technical data:
Width: 150 cm
Height: 202 cm
Depth: 60 cm
With 3 base drawers


Two-person wardrobes/linen cupboards*

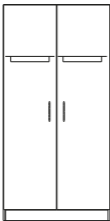


Material no.: 241747
Technical data:
Width: 100 cm
Height: 202 cm
Depth: 60 cm

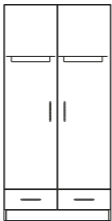

Material no.: 241748
Technical data:
Width: 100 cm
Height: 202 cm
Depth: 60 cm
With 2 base drawers

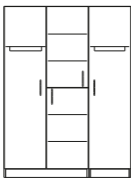

Material no.: 243885
Technical data:
Width: 150 cm
Height: 202 cm
Depth: 60 cm

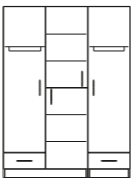

Material no.: 241750
Technical data:
Width: 150 cm
Height: 202 cm
Depth: 60 cm
With 2 base drawers

Extension elements



Material no.: 241738
Technical data:
Width: 50 cm
Height: 202 cm
Depth: 60 cm



Material no.: 241739
Technical data:
Width: 50 cm
Height: 202 cm
Depth: 60 cm
With 1 base drawer

Tops




Material no.: 241754
Technical data:
Width: 100 cm
Height: 42 cm
Depth: 60 cm




Material no.: 241755
Technical data:
Width: 150 cm
Height: 42 cm
Depth: 60 cm

Accessories

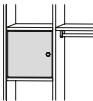

Wooden safe door
Plastic safe box
(without illustration)

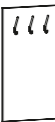

Coat rack
Material no.: 197373
Technical data:
Width: 50 cm
Height: 103 cm



Mirror
Material no.: 197372
Technical data:
Width: 50 cm
Height: 103 cm



Combination Mirror/coat rack
Material no.: 197374
Technical data:
Width: 50 cm
Height: 150 cm
Depth: 22 cm

Decors, tables, coat rack, mirror



Natural Beech R 24010 VV



Royal Maple R 27001 VV



Real Cherry R 42002 VV



Havanna Cherry R 42006 VV



Bella Noce Choco R 42006 VV



Lindberg Oak R 20021 VV

EXAMPLES OF MODELS



Pedestal table, round,
available in a range of diameters



Dining table, rectangular,
available in various sizes

Additional features for tables:
felt, plastic or metal gliders



Coat rack
(W: 50 cm, H: 103 cm)



Mirror
(W: 50 cm, H: 103 cm)

Delivery times
available on
request

*If sideboards or cupboards are combined, the adjoining top surface of the sideboard ends flush with the sideboard body.



Sitting up and taking part in everyday life

Chairs

ADVANTAGES

An attractive range for residences
Burmeier's chairs are comfortable, elegant and practical. They are just as suitable for residents' rooms as for communal areas and waiting rooms in residences. Ergonomic seat heights and armrests offer secure support to people in need of care. They are able take part in daily activities at eye level and maintain their mobility.

Space-saving stacking chairs
Burmeier's stacking chairs can be stacked together quickly to save space. This is practical for floor cleaning as well as when reorganising the room for joint events.

Further information:



The Residenz chair range



Stacking chair
(W: 46 cm, H: 81 cm, D: 56 cm, seat height: 45 cm)



Stacking armchair
(W: 56 cm, H: 81 cm, D: 56 cm, seat height: 45 cm)



Mid-back armchair
(W: 60 cm, H: 91 cm, D: 58 cm, seat height: 45 cm)



High back armchair
(W: 60 cm, H: 104 cm, D: 56 cm, seat height: 45 cm)



Armchair
(W: 57 cm, H: 94 cm, D: 56 cm, seat height: 45 cm)



High back armchair
(W: 57 cm, H: 110 cm, D: 59 cm, seat height: 45 cm)

Additional equipment for seating:
felt, plastic or metal gliders

Minimum order
quantity: 4 chairs;
delivery times
available on request

Burmeier GmbH & Co. KG

Pivitsheider Straße 270
32791 Lage
Phone +49 (0) 5232 9841 - 0
Fax +49 (0) 5232 9841 - 41
info@burmeier.com
www.burmeier.com

A Stieglmeyer-Group company

Hauptsitz / Headquarters

Deutschland / Germany

Stieglmeyer GmbH & Co. KG

Ackerstraße 42
32051 Herford
Phone +49 (0) 5221 185 - 0
Fax +49 (0) 5221 185 - 252
info@stieglmeyer.com
www.stieglmeyer.com

France

Stieglmeyer SAS

8 Allée des Chevreuils
69380 Lissieu
Phone +33 (0) 4 7238 - 5336
Fax +33 (0) 4 7238 - 5481
contact@stieglmeyer.fr
www.stieglmeyer.fr

South Africa

Stieglmeyer Africa (PTY) Ltd

The Estuaries · Century City
Cape Town, 7441
Phone +27 (0) 21 555 - 1415
Fax +27 (0) 21 555 - 0639
info.southafrica@stieglmeyer.com
www.stieglmeyer.co.za

Nederland / Netherlands

Stieglmeyer BV

Peppelenbos 7
6662 WB Elst
Phone +31 (0) 345 581 - 881
info@stieglmeyer.nl
www.stieglmeyer.nl

Suomi / Finland

Stieglmeyer Oy

Ahjokatu 4 A
15800 Lahti
Phone +358 (0) 3 7536 - 320
Fax +358 (0) 3 7536 - 330
stieglmeyer@stieglmeyer.fi
www.stieglmeyer.fi

Polska / Poland

Stieglmeyer Sp. z o.o.

86-212 Stolno
Grubno 63
Phone +48 56 6771 - 400
Fax +48 56 6771 - 410
info.stolno@stieglmeyer.pl
www.stieglmeyer.pl

